


### Henley Treatment and Examination Chair

The Australian designed and made treatment chair is a versatile addition to surgeries, first aid rooms or anywhere there is a benefit in having a patient sitting for treatment. While this chair will also extend to work as a horizontal examination couch, it will easily convert into a comfortable and easily accessible chair. This is achieved through an integrated seat and backrest where, as the backrest lifts, the seat section progressively reclines to -11°. The leg rest has a built-in extension which allows for taller patients and the arm-rests swing back for chair access or to provide easy access when the patient is lying down. Optional wider arms make this chair suitable for blood taking / phlebotomy. The retractable wheel system makes moving the couch easy without compromising its stability in use.



#### Specifications

Product specifications	
Product weight capacity	200 / 350kg
Product length	1930mm
Product width	1000, 1100 , 1200mm
Product height	450 to 950mm
Product weight	85kg
Material	Powder coated steel / foam
Frame colour	Light grey
Vinyl colour	Std blue – others available
Warranty	5 years parts only

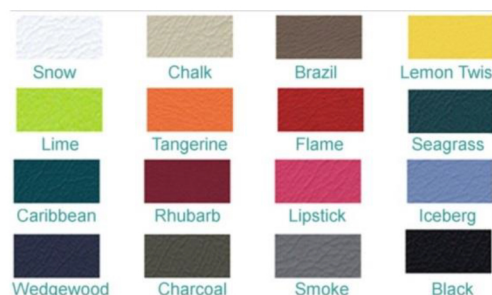

Vinyl cover options