



Active Mobility

Hygiene



Baths



Change tables



Shower trolley



Shower System

Amber

High / Low Change Table

The Reval High/Low Change Table is a wall mounted, height adjustable change table. This change table enables users with poor sitting balance to be ergonomically positioned and changed. The change table is attractive, user-friendly, easy-to-maintain, and has an extensive height adjustment range of 600mm.

The change table is safe for both the patient and the user. The adjustable back rest and fold down safety rail are standard features.

The wall mounted unit is height adjustable for multi-functional use. By positioning this highly versatile change table to the correct height, transfers can be easily made from other aids such as beds, stretchers, shower seats and wheelchairs.

Specifications



Max weight capacity	150kg
Weight of the shower stretcher	85kg
Length of the shower stretcher	1900mm
Wall frame	Steel, electrolytically galvanised
Height adjustable profiles	Stainless steel powder coated
Suspension point for stretcher	Stainless steel
Lifting capacity	135kg
Height adjustment	600mm
Speed	Minimum of 14m0m/s at max load
Power consumption	120W
Drive operates at low tension	24VDC
Protection factor	IP 66

Leichhardt

Height adjustable, junior change table

The Leichhardt Junior Change Table is designed and made in Australia. Ideal for small and growing children, as the cot ends can be used to extend the top from 1250mm to 1650mm long. The minimum height of just 450mm makes access and patient transfers easy and the maximum height of 950mm is high enough for taller carers to work comfortably.

Leichhardt Change Table

Max weight capacity	170kg
Product height	350mm - 950mm (of support surface)
Product length	1250mm
Product width	710mm
Fabric	Medical grade heavy duty vinyl
Foam	38mm high density/12mm soft
Frame	Steel plate & tube in light grey
Wheels	125mm braked
Height adjustment	450 - 950mm



Kensington

Adult Change Table



The **Kensington Change Table** is designed and made in Australia. With a full sized table top, it is long enough and wide enough for most patients.

The minimum height of just 450mm, makes access for frail or disabled patients easy, and the maximum height of 950mm is high enough for tallest carers to work comfortably. Accessory cot sides and cot ends can also be fitted. The four individually braked wheel system makes moving the couch easy without compromising stability in use and provides 155mm of clearance below the frame for patient hoist access. One wheel has a directional locking brake to help control the table when moving.

Features	
Couch length	1930mm
Couch width	710mm
Frame material	Heavy steel plate and tube
Frame colour	Light grey
Height adjustment (lift mechanism)	450 - 950mm
Max weight capacity	200kg
Product weight	70kg

Nivano

Height Adjustable Change Table

The Nivano Height Adjustable Change Table is an easy-to-operate multi-purpose height-adjustable shower/changing table.

Available in lengths from 1200 to 1900mm, the Nivano wall mounted changing bench has a travel range of 700mm, and can easily lift a client from near floor level to a convenient height for washing and changing.

The Nivano folds away very neatly when not in use, protruding just 255mm from the wall and features a battery back-up for safe lowering in the event of a power failure. Available with a range of options including tilting head/back rests, padded tops, guards and protective bumpers; allowing the Nivano table to be tailored just for you.



Specifications

Product lengths	from 1200mm (Paediatric) - 1900mm (Adult)
Max weight capacity	200kg
Height adjustment	700mm
Folded width	220mm

Product features

Wall mounted foldable unit | Durable non-rip high tension vinyl material | Suitable for use in shower areas and wet rooms | Protrudes just 255mm from wall when stowed | Perforated surface for easy drainage | Optional Tilting Back Rest | Height adjustment of 700mm

Nursing Bench

Adjustable Height Nursing Bench

1800mm Powered Nursing Bench with wired hand control, 700 mm height adjustable.

Including safety rail and water collection tray with integrated outlet.

Materials: Coated steel, mesh cover, polystyrene.

The flexibility of the nursing bench provides optimal working conditions for the carer, as well as a high level of comfort and safety for the user.

The integrated gas cylinders enable the bench to be folded out of the way when not in use to maximise floor space in front of the unit.



Specifications

- Remote controlled, low voltage (24V) for safe operation in a wet area
- 700mm of height adjustability (300 - 1000mm)
- No cold surfaces - covered with polyurethane foam
- Curved for added comfort and safety
- Manually adjustable back support
- Water collection tray with integral outlets prevents splashing on carer
- Integrated gas cylinder enables bench to be easily folded against the wall to maximise floor space
- Weight loading capacity 200kg
- Fixed height option available
- Includes safety rail

“The Nursing bench is perfect for aquatic centres, hospitals, nursing homes, retirement villages and schools.”

MONTAR

Fold-up Change Table

Montar Fold-Up Change Table is a wall mounted examination table perfect for small spaces. When folded up, it needs very little space and leaves your examination room free for other uses. It has an adjustable height range of 600mm, and the upholstery is industrial grade vinyl, so spills and splashes are not a problem. This is the perfect examination table for patients with disabilities. The Montar change table also offers remote carer to patient assistance by means of a cordless handset.



Child: Change table / with 3/4 length support surface

Max weight capacity	170kg
Product height	350mm - 950mm (of support surface)
Product length	1420mm (of support surface)
Product width	720mm (of support surface)

Adult: Change table / with full length support surface

Max weight capacity	170kg
Product height	350mm - 950mm (of support surface)
Product length	1900mm (of support surface)
Product width	720mm (of support surface)

SPECIFICATIONS

Electric Hi/Low
Wall Mounted
Standard Blue Padded Vinyl Top
Anti-Bacterial Polyester Powder Coating
Battery Backup
Infra red handset
Drop Safety Rail
SWL 175kg
Low position - 350mm
High position - 950mm
Overall width - closed 350mm Open 1030mm
Overall length - 1880mm
Sleeping surface - 1880 x 700mm



Crystal

Electric or hydraulic shower trolley

The Reval Crystal Shower Trolley is a height-adjustable shower trolley designed to aid with hygiene procedures. The trolley has a mattress, that covers the side rails and a pillow for the comfort of the user. The metal side rails, lock into position for safety and can be easily folded away.

Water can be drained from the trolley with a two metre flexible waste pipe.



Specifications

Product dimensions	Vary with model
Product weight	60kg (hydraulic) and 80kg (electric)
Max weight capacity	160kg (hydraulic) or 150kg (electric)
Sizes	Junior, Adult and Bariatric
Additional information	Hydraulic or electric height adjustment option

Iris

Variable Height Profiling Shower Table

The Iris has been designed to facilitate showering, bathing and hygiene routines for patients; with a wide range of postural and cognitive challenges. Iris also has an ergonomic benefit for care providers, due to its variable seat and platform height. Care workers can maximise their comfort and safety by working at a height which suits them best, ensuring their workplace wellbeing.

Iris can also be used as a mobile transfer system to transport patients and service users in a supine or seated position.



Features

- Multi-position posture management
- Variable height mattress platform
- Comfortable profiling mattress
- Waterproof mattress with high sides
- Lockable side security rails
- Comfort pillow and body belt
- Retracting Stabilising chassis



Waterproof mattress with high side



Iris shown in supine position



Lockable security side rails

Lambda

Mobile Shower Trolley

Lambda shower trolley: Designed to help meet your client needs in the bathroom, whilst also providing a good working environment for the carer, by helping to eliminate heavy lifting.

Lambda is extremely easy to use since all functions are combined in the sturdy hand control. The hi-lo action uses a battery powered motor, giving both a large lifting range and a very smooth action. A separate motor allows the showering surfaces to tilt (Anti Trendelenburg) up to 10 degrees thus ensuring client comfort and adequate water drainage.

- Large and easy to manoeuvre wheels
- Central braking, straight steering and free position
- Adjustable, lockable side and end rails
- End rails have 4 lockable positions



Specifications

Width across side rails	780mm
Length across end rails	2060mm
Width across wheels and centre/centre	660mm and 520mm
Length/ width of base plate	1870/600mm
Product weight	115kg
Max weight capacity	175kg
Lifting interval	560-1030mm (to top of base plate)
Sideways deflection	approx. 450mm in top position
Anti-Trendelenburg angle	approx. 10 degrees
Wheels	150mm
Power supply	24V Linak Jumbo system
Material	Powder coated steel



Onyx Bath

A truly relaxing bathing experience

The Onyx bathing system has been specifically designed for medium to very high dependency needs and is available in two tub sizes. The longer option generously accommodates larger patients, while the compact shorter option is designed to maximise space.

Onyx features a silent high-low function, which allows caregivers to work in a safe position when bathing high dependency patients. HydroSpa, AromaSpa, ColourSpa and AudioSpa sensory options are also available reducing stress, increasing healing times and promoting comfort and relaxation for those patients who need it most.

Specifications

Max water volume	200l (long bath), 180l (short bath)
Adjustable height	620 - 1020mm
Filling time	approx. 5 minutes to 3 bars
Heater	300W
Empty product weight	120kg (long bath), 110kg (short bath)



Transcare Bath

Push-button posture management

The Reval Transcare is a highly versatile bathing system, which is ideal for patients with reduced posture control. The vertical-raise side door allows easy patient access, and the high-low function prevents risk of back strain to the caregiver. The Transcare system also has a unique tilting function, this allows patients with challenging conditions to semi-recline while fully immersed neck deep in water.

HydroSpa, AromaSpa, ColourSpa and AudioSpa sensory treatment options are available with the Transcare, enabling patients to benefit from the therapeutic benefits of massage, colour, music and aromatherapy while enjoying a warm relaxing bath.

Specifications

Vertical travel	150mm
Water reserve	25 litres
Max water volume	90 litres
Max recline angle	22.5 degree
Max weight capacity	200kg
Product weight	140kg



Reinforced fibreglass bath white colour | Handgrip |
Door lock in stainless steel | Electric supply: 240VAc 50Hz |
4 function hand control (up/down and recline/upright) | Back-up battery power

Cocoon

Variable Height Bathing System

With the new **Cocoon Therapy Bath** it is possible to transform a patient's bath into an experience of well-being and relaxation, as well as providing bodily hygiene care.

Select from the user friendly R'GO tub shape with non-slip functionality or the spacious keyhole tub shape.

The tub is available in 4 different colours. Other options include choice of handsfree digital or thermostatic mixer with therapeutic features to enhance the sensory experience – hydrotherapy, aromatherapy, chromotherapy & music therapy.



Features

- Reinforced fibreglass tub – white interior
- Stainless steel handgrips
- Adjustable height
- "anti-vandal" LED safety switch
- Manual waste control
- Overflow outlet
- Digital temperature display
- Compatible with chair and stretcher lifts
- 40mm draining plug with siphon

Caleche

Side access, semi-assistive baths

The **Caleche** is the ideal side transfer bath for either independent users or those requiring assisted bathing. They are suitable for use in hospitals, sheltered accommodation and nursing and residential care facilities and sheltered accommodation.

The **Caleche** and **Caleche Compact** are both stylish and functional, enabling walk in access, and side transfer from a wheelchair or commode chair, as well as from mobile or ceiling track hoists. Their compact designs enable installation in the smallest of spaces, and their large, internal dimensions provide ample room for users up to 200kg. The large, lockable, side access door provides added security. HydroSpa is available as an option.



Specifications

Max water volume	205 litres
Maximum water level	1000mm
Product weight	85kg
Max weight capacity	2300kg

Caprice

Bathing in comfort | Bathing with ease | Effortless bathing

Fixed Height

The highly affordable fixed height Caprice is ideal for users who have reasonable mobility and are able to step in and out of the bath independently or with some assistance. If users need to be hoisted, there is plenty of room under the bath to allow for a mobile hoist.

Variable Height

The stylish variable height Caprice takes the aches and pains out of daily bathing routines. This model features a push-button high/low adjustment, for a safer, more comfortable nursing height.

The height adjustment also enables the bath to be lowered to allow those with reasonable balance and mobility to step in and out. If users need to be hoisted, there is plenty of room under the bath to allow for a mobile hoist.

Variable Height with IntegralHoist and Transfer System

The advanced Caprice, with variable height adjustment and an integrated hoist and transfer system, features the latest bathing technology. Designed to ease the stress and strain of the daily transfer and bathing routine, the user can easily be transferred in and out of the bath, and the tub height adjusted at the press of a button. This makes light and dignified work of transfers, reducing the risk of back ache and pain for the care worker. The low profiling integral hoist saves space in the bathroom, allowing more room to work and eliminating the serious trip hazard of a mobile hoist.

The transfer system and hoist are made from non-corrosive stainless steel and can also be used as a wet room shower chair or over a toilet commode.



Bath

- Reinforced fibreglass bath and frame
- Height adjustable via 2 x 24V electric actuators
- IP67 rated lite touch control
- Integral shower system

Lifter

- Steel chassis
- Polyurethane foam seat and backrest
- Stainless steel seating system and transfer chassis
- Four pivoting brake castors Ø 125mm
- Polyurethane foam folding footrests

Rubis

Mobile Shower System



The new **Rubis mobile shower system** has been designed to shower patients who cannot access the bathroom independently or with assistance.

The Rubis is compact and highly manoeuvrable, it is also supplied with a purpose designed, double-width drying towel and waterproof shower sheet. All components are fully integrated into the appliance making it safe and user friendly.

Features

- Fresh water tank in stainless steel with a capacity of 23 litres
- Waste water tank in stainless steel with a capacity of 23 litres
- Temperature controller with double display
- Red warning light if water is too hot
- Handheld shower with flow adjustable trigger
- White anti-corrosion coated steel frame
- Six directional wheels, two of which has brakes



Control panel



Fresh water tank



Mobile shower system

Alto Commode

Care Chair

The Alto CareChair is the only integrated care chair where all functions are electrically operated. Lifting, seat tilt, central braking and even directional lock are all accessible via the hand control. Despite an impressive maximum user weight of 150kg, the Alto CareChair is incredibly compact and fits into the smallest of bathrooms. With an electrical lifting range of 470mm it is possible to transfer a patient directly from the bed without a patient lifting hoist.

Specifications

Product height	1070 - 1535mm
Product width	570mm
Product length	980 - 1290mm
Seat height	520 - 985mm
Seat width	450mm
Seat depth	430mm
Turning diameter	1040mm
Product weight	44.5kg
Lifting machinery	Electrical actuator
Battery capacity	25 - 30 lifting cycles
Max weight capacity	159kg



Hera

Tilt-in-space bathing & toileting chair

The Hera bath and toilet chair, with lift up arm rests, has been developed with the carer and user in mind. The tilt function is easily activated using the cable under the handle. The chair can be placed in any desirable position using the gas spring i.e. in a horizontal or tilt position.

A unique feature of this chair is the built-in height adjustment, which can be carried out without using any tools. The seat is made of black PU foam, and is available with or without the front opening. The back support, also made of PU foam, is height adjustable for better comfort.

Specifications

Max weight capacity	140kg
Product length	710mm
Product width	580mm
Product height	500 to 750mm
Product weight	35kg
Material	Steel / PVC
Colour	Grey / Black



Rise Chair

Mobile Shower System



The Rise IIS-130 integrated care chair is incredibly compact and able to fit into the smallest of bathrooms. With an electrical lifting range of 420mm, carers are able to transfer a patient directly from the bed without using a patient lifting hoist. This lifting function also allows the patient to use the basin and mirror effectively, thus creating a more dignified care routine. The electrical seat tilt (+6°/-31°) and multi-adjustable soft backrest offers increased user comfort and stability.

Specifications

Product height	1080 - 1480mm
Product width	570mm
Product length	880 - 1050mm
Seat height	535 - 950mm
Seat width	450mm
Seat depth	430mm
Turning diameter	980mm
Product weight	32kg
Lifting machinery	Electrical actuator
Battery capacity	25 - 30 lifting cycles
Max weight capacity	130kg

Soflex

Shower and Toilet Chair



The Soflex Shower and Toilet Chair is an ergonomic solution for the care and hygiene routines of patients. It offers optimum comfort, for both users and carers. The Soflex can also be used as a commode chair, used with a bowl or pan, or directly over the toilet.

Available with hydraulic or electric controls, the Soflex is height adjustable to enable carers to work at a height that suits them best. The tilt-in-space function allows for use by patients with limited upper body stability.

Features

- Head and footrest that can be adjusted to match all heights and body types
- The specially shaped and comfortable base facilitates washing, easily fits a WC pan or provides direct access to the toilet seat
- Takes up little floor space and features four swivel wheels with brakes
- Available in hydraulic or electric version

XXL Care Lift

Bariatric Commode

The XXL Care Lift is a completely new concept in the area of bariatric care. With a 600 mm lifting interval, as well as electrical tilt-in-space and elevating leg support; the XXL Care Lift offers ideal working conditions for caregivers during toileting, showering and general care routines. For patients weighing up to 325kg, this product offers optional headrests in two sizes, a multi-adjustable soft backrest and also leg supports that can be easily folded away, for improving access for ambulant patients.

Specifications

Product weight	325kg
Product width	1350mm
Product depth in retracted position	780mm
Product dept in extended position	1230mm
Seat height	500 - 1100mm
Seat width	770mm
Seat depth	540mm
Product weight	120kg
Turning diameter	1300mm
Lifting mechanism	24 V electrical with 10 - 15 lifting cycles



Topaz

Bath Stretcher System



Our Topaz Bathing Transfer Stretcher makes bathing efficient and comfortable. It features three wide anti-shear mattress sections, which fit neatly into variable height baths, as well as domestic baths, subject to installation. The mattresses are easily removed for cleaning and replacement, the safety armrests fold down, the leg and backrest sections are fully adjustable, and the anti-corrosion steel frame has been specifically designed for full immersion.

The Topaz stretcher provides effective posture management for patients with more challenging needs, and is easily manoeuvrable with the swivel wheels and low wheelbase, allowing precise positioning in to the bath.

Available in powered or manual options

Features

- 24V electric actuator
- Electric remote control, emergency stop button, emergency lowering button
- Range 80 - 130kg
- Hydraulic actuator (Manual lifting system only)
- Foot operated adjustment

Topaz

Bath Hoist System



The Topaz bathing Hoist makes bathing routines easy. Compact in design, and light and easy to manoeuvre, the Topaz is ideal for use with height adjustable baths. The low wheelbase, and four swivel wheels and brakes, enable precise positioning underneath the bath.

Patients can be immersed into and out of the bath effortlessly, smoothly and silently at the press of a button. The Topaz also features a steel structure that is corrosion resistant, and a safety belt, providing added security. The armrests and padded backrest fold away and the comfortable one-piece moulded seat is removable for cleaning, ensuring maximum hygiene during the healthcare routine.

Available in powered or manual options

Features

- 24V electric actuator
- Electric remote control, emergency stop button, emergency lowering button
- Range 80 - 130kg
- Hydraulic actuator (Manual lifting system only)
- Foot operated adjustment

Poolside Hoist

Pool lift

The Reval Pool Lift is a deluxe pool access system allowing autonomous access to the pool thanks to the rotation system and the electrical hi-lo, which is piloted by a pneumatic remote control. The revel is fixed to the ground, meaning this pool hoist allows the safe transfer of handicapped people into chairs or laying down on a stretcher.

Specifications

Product weight	130kg
Max weight capacity	150kg
Max immersion of the patient	63cm
<ul style="list-style-type: none">• 2 batteries• 1 recharger, recharging time: 6 hours• 1 green luminous indicator "ON" indicating the pool lift is under tension• 1 battery level indicator• Submersible pneumatic control, Emergency stop button,• Anti-entrapment device, Emergency rise and lower buttons,• Rotation cut out device for anti-impact and entrapment	



Penguin Hoist

Pool lift

Designed and manufactured in Australia, the Penguin Pool & Spa Access Hoist allows a person living with any disability to have access to both above and below ground Spas, Swimming & Hydrotherapy pools, without any manual lifting. It is suitable for use in public pools, leisure centres, private homes, nursing homes and schools. The unit can be mobile, or fixed permanently at the edge of the Pool or Spa. It is easily installed on the pool deck by a handy person, or alternatively we can arrange installation for you. All pools are different, at AMS we will work with you to ensure the Penguin Pool and Spa Hoist meets your pools needs.

Specifications

Max weight capacity	150kg and 200kg
Batteries	Linak- 24V re-chargeable battery
Rotation	Manual - 360degrees
Product weight	60kg
Floor mounting plate	450 diameter
Construction	Stainless steel, powder coated to colour of choice
Hand control	Push button
Installation	Anchored to pool deck with Stainless Steel



Guardian

Bath Stretcher System



The Guardian Pool Hoist is easily operated by the attendant, with a hydraulic-powered lifting mechanism, which allows for recreational and therapeutic swimming. This makes the experience more enjoyable for both patient and attendant. The stainless steel painted finish is designed to resist the humid and corrosive swimming pool environment, and the heavy duty construction supports up to 182kgs.

Features

Product weight	80kg
Max weight capacity	182kg
Mounting style	Cement-in sleeve/ Surface mount
Mast Height cement-in sleeve	1295mm
Mast Height surface mount	1524mm
Range of boom cement-in sleeve	620 - 2190mm
Range of boom surface mount	845 - 2420mm
Boom length	1150mm
Finish	Grey stainless steel paint

Hydraulic power | Extension arm

Reval Butterfly

Hydrotherapy Pool

The Reval Butterfly Pool provides relaxation for up to four people. It is a cost effective solution when upgrading your facility because it is a surface mounted pool that can be installed virtually anywhere. Designed to stand-alone, it allows the caregiver 360 degree access, and comes with a stainless steel surrounding handrail for safety. A range of accessories and access solutions complement the Reval Butterfly Pool.

Specifications

Water capacity	1800 litres
Overall length	2480mm
Width of basin	2000mm
Height	900mm

- 4 handrails in stainless steel
- Blue colour finish, in reinforced acrylic
- Metallic frame in stainless steel - equipped with 8 adjustable feet for height adjustment
- Thermostatic mixer and cleaning shower



Reval Tonic

Swim and gym in one

The Reval Tonic Hydrotherapy Pool is at the cutting edge of rehabilitation hydrotherapy. Available in several sizes, it has a full-length level access platform that lowers into the pool to a suitable depth.

The pool incorporates a dynamic hydro flow system, and a water current control function which allows the patient to exercise in moving water; helping build movement, muscle tone and strength.

Specifications

- The Tonic pool affords up to three users at any one time
- Its adjustable walk-on platform provides safe and easy level access for users
- The platform can be lowered to a desired depth to meet varying patient's needs at the press of a button
- The Tonics filtration system is safe and effective. The filtration cycle is hourly. The entire pool volume of water is filtered and cleaned to ensure maximum hygiene
- Tonic is available complementary range of submersible accessories: Guide and hand rails, Aqua wheeled chairs, stand-up aqua trolleys, etc...
- Tonic is available with hydro massage system designed to further promote relaxation and to build muscle tone.
- Easy access



Balnaeo Pool

Therapy at its most flexible



The Balnaeo hydrotherapy pool is a made to measure pool system. Clients can specify their requirements and Reval will meet their functional and architectural needs.

Supplied with a specifically designed plant room, the Balnaeo assures properly calibrated filtration and chlorination, making routine pool maintenance easy.

The pool is available in variable depths, and with full or sectioned inbuilt hydrotherapy systems. Several in-pool apparatus are available to support a range of therapy and exercise regimes.

POOL ACCESSORIES

In-pool exercise & parallel bars

Ideal for underwater physiotherapy and exercise routines, our in-pool exercise bars can be used with many techniques and exercises that strengthen and build muscle tone. Made from high quality stainless steel, the exercise bars are perfectly designed for rehabilitation work, and are easily repositioned along the poolside for exercising at different water depths.

Alternatively, our in-pool parallel bars can accommodate multiple users. With an option of one or two walking lanes, the parallel bars are fixed to the pool floor, providing submerged support for load bearing exercises and can be manufactured to specific requirements



In-pool exercise chair

Our in-pool exercise chair is a useful underwater aid. The submerged seat can be easily repositioned along the pool's side rail, easily allowing those with posture control challenges to follow exercise routines without the need to stand or swim.



Hydro-Massage Lance

For direct muscle stimulation and improved circulation, our submergible hydro-massage lance assists underwater physiotherapy routines. The temperature is set to a comfortable 36 degree, while the water jet intensity and pattern can be adjusted when working on particular parts of the anatomy.



In-pool steps

Built to any width or height (in accordance with regulations) our poolside and in-pool steps are made from high grade stainless steel, offer a safety handrail and easy clean, removable, polyurethane walking surfaces.



Aqua Chair & Stretcher

Operated by a user-friendly waterproof hand control, the poolside hoist smoothly and quietly lifts the aqua stretcher or chair, leaving the chassis and wheels at the poolside. When submerged the stretcher platform tilts, easing the patient into the water weightlessly and with little effort. The pool stretcher may also be used in conjunction with our level access platform.

Both have easy clean surfaces, the aqua chair has a comfort moulded seat and the aqua stretcher has a padded mattress and fold-down safety sides.





Tiltlock Folding Grab Rail

The Tiltlock is a wall mounted, double arm, folding grab rail that automatically locks when folded down to the horizontal position, and locks vertically out of the way when not required. It is designed to be positioned on either or both sides of a toilet to assist with transfers. It is available in one size and four finishes.

An optional toilet roll holder and nurse call mount block is available (nurse call switch is not included).


Length of rail: 850mm



Grab Rails

A full range of standard sizing and custom made grab rails are available in powder-coated stainless, brushed stainless and copper finishes.

Please enquire for further information and specifications.

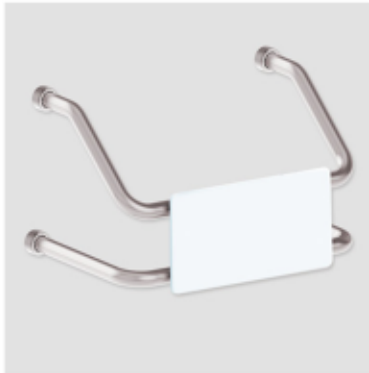

Safe-Assist Backrest

The Backrest, a support device for toilets, bidets and the like, has been introduced into the Australian Standard for Access and Mobility to comply with the changing requirements for sanitary facilities.

The Backrest provides a support that improves the comfort of using a toilet, while simultaneously providing an acceptable level of safety for disabled, aged and infirmed.



Fold-away shower seats

When safety in the shower is paramount, shower seats are a solution; with rounded edges and seats that are foldable, self-draining and non-slip. Manufactured from 304 grade stainless steel, a white, high density polyethylene seat top and also fitted with anti-drop tension hinges that ensure stability up to 160kg.

Shower Chair

- Height adjustable
- Ergonomic angled arms assist transfers
- Drainage slots drain water away from seating area
- Non-slip rubber tips
- Electrophoretic surface treatment to prevent rusting
- Powder coated steel frame

Specifications

Seat Width	410mm
Seat Height	400-550mm
Max Base Width	560mm
Max Base Depth	580mm
Max User Weight	125kg
Width Between Arms	460mm



Shower Stool

- Suits smaller shower areas
- Height adjustable
- Integrated arms assist transfers
- Durable plastic moulded seat
- Non-slip rubber tips
- Electrophoretic surface treatment to prevent rusting
- Powder coated steel frame

Specifications

Width Between Arms	460mm
Seat Depth	310mm
Seat Height	510-660mm
Max Base Width	570mm
Max Base Depth	410mm
Max User Weight	125kg



Overtoilet Aid - Aluminium

- Provides a raised seat and support arms to assist toilet transfers
- Ideal for people with hip or knee injuries, or those who require greater seat height
- Ergonomic angled arms assist transfers
- Height adjustable
- Removable toilet seat and splashguard
- Non-slip rubber tips
- Powder coated aluminium frame is lightweight, rust proof and ideal for wet environments
- Tested and conforms to Australian Standards AS/NZS 3973: 2009

Specifications

Seat Height	470-620mm
Width Between Arms	445mm
Arm Height	230mm
Max Base Width	515mm
Max Base Depth	620mm
Max User Weight	125kg



Disposable Slings



The single patient use sling can be used as a general purpose sling for short stay at hospital or residency visits. Featuring stiffened head support, easy to apply and all supplied with padded legs and additional tailoring to commode aperture, as standard. Used specifically in acute markets to limit cross infection.

Comes in a **box of 10**
Available in **loops** or **clips**
SWL 200kg

Care Sling



The care sling is excellent for transfer onto a toilet, as it provides an open area from the middle of the back to the middle of the thigh. It is recommended for general toilet transfers.

The Care Sling should only be used if the patient has sufficient trunk and upper limb control to maintain balance. The top of the sling must pass below the armpits to maximise stability. Excessive pressure under the armpits should be avoided.

Head Support Care Sling



Head support care sling is excellent for transfer onto a toilet, as it provides an open area from the middle of the back to the middle of the thigh. It is recommended for general toilet transfers.

The Head Support Care Sling is ideal for use if the patient has limited trunk, upper limb and head control. The sling provides greater stability and support for upper trunk and head areas.



Active Mobility

Active Mobility Systems

Unit 8, 110-120 Silverwater Rd., SILVERWATER NSW 2128

P: 02 9649 2111 F: 02 9649 8506

W: www.activemobility.com.au E: sales@activemobility.com.au