

# MOBILITY SCOOTERS





Do you need a lightweight scooter that can fit easily into the boot of a car?

Y

N

Do you need to travel longer distances?

Y

N

**LEO™**  
THE COMPACT SOLUTION



**COLIBRI™**  
THE PORTABLE SOLUTION



Do you need a nimble, mid-size scooter for indoors and outdoors?

Y

N

**PEGASUS™**  
THE MULTIPURPOSE SOLUTION



Will you need support of up to 220kgs?

Y

N

**COMET HD™**  
THE HEAVY DUTY SOLUTION



Do you need a large powerful scooter?

Y

**COMET ALPINE™**  
THE ALL-TERRAIN SOLUTION



# WHICH SCOOTER IS RIGHT FOR YOU?



*Yes, you can.®*



# FREE 24 Hour Roadside assistance program



## Who Are 24/7 Roadservices?

24/7 Roadservices was established 10 years ago as a Roadside Assistance package for new and used car warranties. Following the success of Roadside Assistance as a backup for warranty, they decided to launch the program to the general public and the wholesale automotive industry.

## Frequently Asked Questions

### What does the 24 hour Roadside Assistance program provide Invacare customers?

.....

This program is not a repair service, it is simply a transport service for you and your mobility scooter. If something goes wrong with your mobility scooter and you are unable to drive home, just call 24/7 Roadservices on their 1800 number and they will send somebody to pick you and your scooter up, and transport you to an address of your choosing. You might elect to go directly to your Invacare dealer so the scooter may be serviced, or you might decide to go home and arrange to have your scooter serviced another day.

### Nationwide Service Network

.....

Local service and maintenance support is close at hand via Invacare's Nationwide Service Network. You can either contact your local dealer directly or you can call 24/7 Road Services on 1800 008 616 and we will put you in contact with your nearest service agent.

### How do become eligible for the 14 months FREE offer?

.....

Simple. If you purchase any Invacare mobility scooter, you need to ensure your Invacare approved dealer registers for your warranty online at the time of purchase\*. Once your dealer has registered for your warranty, you will be automatically signed up for the program and can call the 1800 number any time you need to!

\*If your dealer fails to register for your warranty at the time of purchase, you can arrange this at a later date however, you will not be registered for the 24 hour roadside assistance program until you are registered for your warranty.

### What if I want to continue with the 24 hour Roadside Assistance program after my 14 months is up and how much will this cost?

.....

No problem! Just phone 24/7 Roadservices to sign up. They will talk you through the process. Annual fees are subject to change over time.

### I live in a remote area. Can 24/7 Roadservices send someone to collect me if I break down?

.....

24/7 Roadservices is a national network of people to call on when you need help. We feel confident that no matter where you live, this service will be available for you.







### Micro, Light, Portable: the *Invacare Colibri*

The new Invacare Colibri is a stylish, simple and colourful micro lightweight scooter, designed for those who enjoy a totally independent lifestyle. Thanks to the unique Invacare LiteLock™ system the scooter can easily be taken apart without the need for tools and fits neatly into any car boot.

- ✓ Fits easily into any car boot. **Great for tight spaces.**
- ✓ Powerful motor. **Travels up to 8km/h.**
- ✓ Full sized swivel seat. **Easy to get in and out of.**
- ✓ Simple dismantling. **For easy transportation.**

**NEW**  
with  
18Amp Battery



#### Effortless dismantling

Don't think - do it! Less than a minute is all that's needed to disassemble the Colibri for transportation.


#### Components weight

Front section	15.6 kg
Drive unit	9.6 kg
Seat	9.6 kg
Battery box 18 Ah	14.1 kg



#### Increased comfort

The soft, padded seat provides comfort during driving. For increased comfort there is also the optional 'easy fit' seat suspension kit.

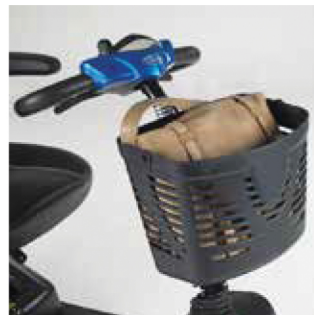

#### Ergonomic Handlebar and Plastic Basket

For greater comfort and prevents fatigue in hands and wrists.



#### Choose from five great colours

## Features

Features	Benefits
<b>LiteLock™</b>	Easy to take apart and assemble
<b>Flat Free Tyres</b>	Prevents scuffing and punctures
<b>Ergonomic Handlebar</b>	Prevents fatigue in hands and wrists
<b>18 Amp Battery</b>	More comfortable, powerful ride

## Technical Specifications



47cm



40cm



370 - 420mm



39cm



51cm



101cm



35kg\* - 49kg^



136kg



2 x 18 Ah



5cm



6°



200W



16cm



8 km/h



220cm

#  
^  
\*

Excluding batteries  
Including batteries  
Peak performance

To add flag, optional  
bracket available.

Optional suspension  
kit available in low,  
medium and high.

## Shroud Colours

	Scooter Colibri white 4w	PKG-1577168
	Scooter Colibri jasper red 4w	PKG-1577169
	Scooter Colibri jade green 4w	PKG-1577170
	Scooter Colibri sapphire blue 4w	PKG-1577171
	Scooter Colibri stone orange 4w	PKG-1577172





### Safe, Solid and Secure...

The Invacare Leo™ is an affordable lightweight scooter ideal for both indoor and outdoor locations.

Planning a short break or weekend away? The Leo conveniently dismantles into five separate parts. The Leo neatly fits into the boot of most cars and can be easily dismantled and re-assembled without tools.

### Standard Configuration

2x 36 Amp batteries  
Off-board 8 Amp charger  
Safety Vest

### Codes

Model PKG-HS-528-BU4 (Sapphire Blue)  
PKG-HS-528-RD2 (Ruby Red)

### Features



2-step disengaging lever prevents scooter from being knocked into freewheel mode



Built in anti-splash guards protect the electronics and transaxle from dirt & water damage



Ergonomic Handlebar and Plastic Basket for greater comfort and prevents fatigue in hands and wrists.

Features	Benefits
<b>Automatic turn-off indicators (after 30 seconds)</b>	Alerts pedestrians to your movements reducing the risk of collisions
<b>Intelligent Lighting Package</b>	If one light bulb blows all other lights continue to function improving visibility
<b>Audio and Visual low battery warning</b>	Alerts you when battery power is low so you are less likely to get a flat battery when out and about
<b>Anti-corrosion chassis</b>	Prevents rusting of the scooter prolonging its life

### Technical Specifications



47cm



41cm



46-56cm



47.5cm



59cm



122cm



59kg<sup>#</sup>- 83kg<sup>^</sup>



136kg

#  
^  
\*

Excluding batteries  
Including batteries  
Peak performance



2 x 36 Ah



131cm



6cm



7cm



10°



240W



36km



8 km/h





### Stable, Trusty and Durable...

With a focus on safety and reliability, The Invacare Pegasus is equipped with multiple features to ensure a safe and easy ride. Combining ergonomics with style, the Pegasus is the obvious choice for an enjoyable, relaxing ride - whether you're travelling near or far.

### Standard Configuration

2x 50 Amp batteries  
Right handbrake  
Safety vest  
Safety flag  
Left and right rear mirrors  
Front basket

### Codes

Model    PKG-1585282 Silver 4W  
            PKG-1585246 Onyx Blue 4W  
            PKG-1585247 Ruby Red 4W

### Features

Features	Benefits
<b>Ergonomically shaped seat with soft padding; recliner lever; and height adjustability</b>	Provides optimum comfort and positioning
<b>Ergonomic handlebars</b>	Prevents fatigue in the hands and arms
<b>Automatic speed reduction on bends</b>	Reduces risk of falls/tipping
<b>Speed reduction button</b>	Provides a simple way to halve current speed at the touch of a button – handy in crowded areas
<b>Brake light with delay</b>	Warns pedestrians that the scooter is stopping reducing the risk of collisions
<b>Handbrake assists the user to stop more suddenly if needed</b>	Reduces risk of a collision in unexpected situations
<b>Tie-down points to prevent the scooter from moving around when being transported</b>	Improves safety during transport*
<b>Anti-corrosion chassis</b>	Prevents rusting of the scooter prolonging its life


















A fully suspended chassis provides a comfortable ride at all times



Built in anti-splash guards protect the electronics and transaxle from dirt & water damage

\*NB: Invacare recommends that, when a scooter is being transported in another vehicle, the user sits in the *vehicle* seat with his/her seatbelt on, NOT on the mobility scooter that is being transported

### Technical Specifications

4W	 51cm	 47cm	 42-49.5cm	 50cm	 62.5cm	 129m	 78kg*/110kg^	 136kg
4W	 2 x 50 Ah	 140cm	 8cm	 9cm	 12°	 240W/600W*	 42km	 10 km/h*

# Excluding batteries  
^ Including batteries  
\* Peak performance





### Powerful, Determined, Adventurous.

New to the Invacare Comet family, the Comet Alpine version offers improved traction and a larger motor, ideal for challenging outdoor uses. The robust power unit and 4-point independent suspension system ensures a continuous comfortable ride even on bumpy ground and can even tackle the most challenging of hills.

### Standard Configuration

Front basket  
2x 75 Amp batteries  
Left handbrake  
2 Rear mirrors  
Safety flag  
Safety Vest  
Rubber mat

### Codes

















Model PKG-1585250 Onyx Blue AL  
PKG-1585251 Ruby Red AL  
PKG-1585306 Silver AL

### Features

Features	Benefits
<b>Larger motor</b>	Enables more power for outdoor use (650W – 1800W)
<b>Full Suspension</b>	A fully suspended chassis provides a comfortable ride even on bumpy ground
<b>Electric seat raiser</b>	Enabling you to electrically raise the height to the preferred height
<b>Off-board charging</b>	Allows you to conveniently charge the batteries away from the scooter
<b>Automatic speed reduction</b>	Automatically reduces your speed around bends preventing you from going too fast when turning and reducing risk of falls

\* NB: Invacare recommends that, when a scooter is being transported in another vehicle, the user sits in the *vehicle* seat with his/her seatbelt on, NOT on the mobility scooter that is being transported

### Technical Specifications

								Compliance: TGA-registered (ARTG# 119338) CE Marked TUV certified – tested to European Standard (DIN EN 12184:1999) * Without batteries
51cm	48cm	44-51.5cm	50cm	66cm	147.5cm	143kg*	160kg	
								
(2)12V 75Amp	163cm	10cm	13.5cm	10°	Cont: 650W Peak: 1800W	50km	10 km/h	

# Invacare® | Comet™ Heavy Duty

**220kg CAPACITY WITH  
HIGH TORQUE MOTOR**



## Well-proportioned, Comfortable, Supporting.

The Invacare Comet HD is part of the Invacare Comet family, renowned for its safe, fast and enjoyable drive. This rugged outdoor vehicle is ideal for bariatrics users weighing up to 220 kg. Packed with safety features, this powerful scooter will take you to your destination in the most direct and secure way.

## Standard Configuration

Front basket  
2x 75 Amp batteries  
Left handbrake  
2 Rear mirrors  
Safety flag  
Safety Vest  
Rubber mat

## Codes

Model PKG-1585252 Ruby Red HD

## Features

Features	Benefits
<b>Larger motor</b>	Enables more power for outdoor use (650W – 1800W)
<b>Full Suspension</b>	A fully suspended chassis provides a comfortable ride even on bumpy ground
<b>Electric seat raiser</b>	Enabling you to electrically raise the height to the preferred height
<b>Off-board charging</b>	Allows you to conveniently charge the batteries away from the scooter
<b>Automatic speed reduction</b>	Automatically reduces your speed around bends preventing you from going too fast when turning and reducing risk of falls

\* NB: Invacare recommends that, when a scooter is being transported in another vehicle, the user sits in the *vehicle* seat with his/her seatbelt on, NOT on the mobility scooter that is being transported

## Technical Specifications



66cm



53cm



44-51.5cm



61cm



66cm



147.5cm



148kg\*



220kg

Compliance: TGA-registered (ARTG# 119338) CE Marked TUV certified – tested to European Standard (DIN EN 12184:1999)

\* Without batteries



(2)12V 75Amp



326cm



10cm



13.5cm



10°



Cont: 650W  
Peak: 1800W



50km



10 km/h



# Invacare® | Superior Safety Features




## Hand brake

Ensures immediate braking if required.



## Automatic speed reduction

This feature automatically reduces your speed around bends.

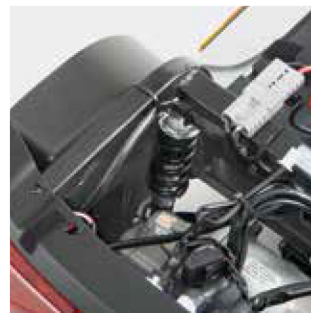

## Full Suspension\*

A fully suspended chassis provides a comfortable drive at all times.

\*Suspension options available for Colibri



## Full lighting system with dual bulb headlights

Dual headlights provide increased visibility. Indicators automatically go off after 20 seconds. Indicators can easily be programmed with or without an audible signal.



## Low battery warning

A visual warning on the battery level display and an audio warning (beep), tells you when the battery needs recharging.



## Integrated brake light

Heavy duty brake lighting warns people around you that the scooter is slowing down.

## Audible reverse alert

An audible signal (beep) lets people know that the scooter is reversing.

## Anti-tip wheels

Prevents the scooter from tipping backwards when climbing steep slopes or negotiating high obstacles.

## Anti-corrosion chassis

Anti-corrosion-treated chassis prolongs the life of the scooter.

## Large wheels

Provides greater feeling of safety when negotiating kerbs and obstacles.



## Two-step disengaging lever

Prevents the scooter from accidentally free-wheeling.

## Finite tiller movement

The tiller can be easily adjusted with a lever to suit your needs.

## Control panel

A clear, intuitive keypad with domed buttons ensures accurate operation. LED and audio signals indicate when each function is activated.

## Ergonomic steering

Prevents fatigue in arms and hands during long trips.



## Built-in anti-splash guards

Protects electronics and transaxle from dirt and water.

## Intelligent lighting system protection

Only the bulb that has blown goes out: all others remain lit (scooter is still visible).

## Dependable electronics

Best-in-class quality wiring and electrical connections optimise reliability.

## Guaranteed reliability

Maintenance and easy serviceability are key design features. Batteries and electronics are easily accessible.



## Accessories and Features

### Adjustable



#### Seat Raiser\*

An electronic mechanism raises the seat to your referred height.

*\* Only available for Comet Alpine*

#### Seat Suspension Module

Adds extra comfort during driving. Kits available for all scooters.

### Safety



#### Off-board charging kit

Simple remove batteries to charge with the off-board charger kit.



#### Seat Belt

Ensures you are kept securely and comfortably in place.



#### Mirror

Moveable mirror can be mounted to the left or right.



#### Flag



#### Hour Counter

allows you to track scooter life by hours operated, perfect for anticipating regular service requirements.



#### Safety Vest



## Storage



### Lockable front box

Secure box for personal valuables  
SP1551638  
*Fits Leo, Pegasus & Comet*



### Rear Seat basket

Perfect for bulky or extra shopping  
1505727  
*Fits Leo, Pegasus & Comet*



### Crutch & Cane holder

Keeps your crutch or cane when on the move.  
1505684  
*Fits Leo, Pegasus & Comet*



### Walker Holder

1505728  
*Fits Leo, Pegasus & Comet*



### Oxygen cylinder holder

Safe storage of oxygen during transit  
1505729  
*Fits Leo, Pegasus & Comet*  
- suitable for B, CS, & CH  
Oxygen Cylinders



### Cup Holder

Stores water or a soft drink  
SP1521594  
*Fits Pegasus & Comet*



### Rear Seat bag

1505731  
*Fits Leo, Pegasus and Comet*



### Heavy-duty Storage cover

Keeps the scooter clean and dry  
1509434 *Fits Leo*  
SP1518240 *Fits Pegasus & Comet*

## Upholstery colour



### Black vinyl

Comet Alpine  
Comet HD  
Colibri, Leo &  
Pegasus

Scooter Comparison	Micro	Small
		
	Colibri	Leo
<b>Model</b>	4 wheel	4 wheel
<b>Colours</b>	White, Red, Green, Blue, Orange	Sapphire Blue, Ruby Red
<b>Overall Length</b>	101 cm	122 cm
<b>Overall Width</b>	50 cm	59 cm
<b>Wheels: Front &amp; Rear</b>	21cm/ 8"solid	26cm/ 10" Pneumatic
<b>Weight without batteries</b>	35 kg	59 kg
<b>Max. Speed</b>	8 km/h	8 km/h
<b>Weight Capacity</b>	136 kg	136 kg
<b>Ground Clearance</b>	3.5 cm	7 cm
<b>Grade Climbable</b>	6°	10°
<b>Kerb Climbing</b>	4.5 cm	6 cm
<b>Turning Radius</b>	220 cm	262 cm
<b>Suspension</b>	Optional	n/a
<b>Seat Type</b>	Ergonomic, mid-back seat	Ergonomic, mid-back seat
<b>Seat Width</b>	46.5 cm	47 cm
<b>Armrest Width</b>	Extension of 5 cm / side	Extension of 5 cm / side
<b>Continuous &amp; Peak Motor Size</b>	Cont. – 200W	Cont. – 240W, Peak – 500W
<b>Standard Battery Size</b>	12V 18H	12V 36AH
<b>Standard Travel Range</b>	Up to 16 km	Up to 38 km
<b>Standard Battery Charger</b>	On-board DC 24V / 5 Amp	DC 24V / 5 Amp
<b>Electronics</b>	DR-50 50A	Dynamic DR-90 90 Amp Controller
<b>Limited Warranty</b>	Scooter 2 years, Batteries & Charger 1 year	Scooter 2 years, Batteries & Charger 1 year




Medium	Large	Heavy Duty
		
<b>Pegasus</b>	<b>Comet Alpine</b>	<b>Comet Heavy Duty</b>
4 wheel	4 wheel	4 wheel
Moonstone Silver, Onyx Blue, Red	Red, Silver, Blue	Red
129 cm	147.5 cm	147.5 cm
62.5 cm	66 cm	66 cm
29cm/ 11" Pneumatic	31cm/13" Pneumatic	31cm/13" Pneumatic
78 kg	143 kg	148 kg
10 km/h	10 km/h	10 km/h
136 kg	160 kg	220 kg
9 cm	13.5 cm	13.5 cm
12°	10°	10°
8 cm	10 cm	10 cm
140 cm	326 cm	326 cm
Front spring; Rear spring & Adjustable	Front & Rear, Spring & Adjustable	Front & Rear, Spring & Adjustable
Extra soft, ergonomic, recline lever, high back with Headrest	Extra soft, ergonomic, recline lever, high back with Headrest	Extra soft, ergonomic, recline lever, high back with Headrest
51 cm	51 cm	66 cm
51-59 cm	51 – 59 cm	51 – 66cm
Cont. – 240W, Peak – 600W	Cont. – 650W, Peak – 1800W	Cont. – 650W, Peak – 1800W
12V 50AH	12V 75AH	12V 75AH
Up to 42 km	Up to 50 km	Up to 50 km
DC 24V / 8 Amp	DC24V/ 8 amp	DC24V/ 8 amp
Rhino 2 120 Amp Controller	Rhino 2 160 Amp Controller	Rhino 2 160 Amp Controller
Scooters 2 Years, Batteries & Charger 1 year	Scooters 2 Years, Batteries & Charger 1 year	Scooters 2 Years, Batteries & Charger 1 year

# Invacare® | Portable Ramps

Invacare Portable Ramps are manufactured with a unique combination of uni-directional glass and carbon fibre materials to create an extremely strong and lightweight ramp suitable for a wide range of uses.

Invacare Portable Ramps have been strength tested to AS/NZ Standards for loads up to 300kg, and require no assembly, so can be deployed for use in seconds!

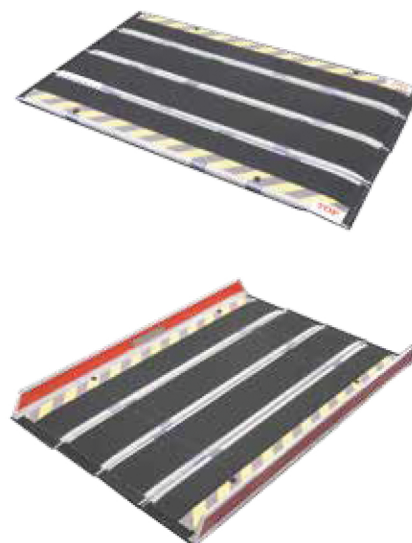

<b>Model/Length (cm)**</b>	<b>Max Height**</b>	<b>Ramp Width</b>	<b>Ramp Weight**</b>
Personal 070	11.5cm	74.5cm	3.9kg
Multipurpose 087	14.5cm	83.5cm	5.8kg
Multipurpose 100	16.5cm	83cm	7.5kg
Senior 090	15.0cm	75cm	6.4kg
Senior 120	19.5cm	75cm	8.0kg
Senior 135	22.0cm	75cm	8.1kg
Senior 150	24.5cm	75cm	10.0kg
Senior 165	13.5cm	75cm	11.4kg
Senior 185	15.5cm	75cm	11.9kg
Senior 200	16.5cm	75cm	12.3kg
Senior 250	20.5cm	75cm	17.3kg
Senior 300	25.0cm	75cm	22.1kg
Senior 350	29.0cm	75cm	24.7kg
EBL 090	15.0cm	75cm	7.4kg
EBL 120	19.5cm	75cm	9.7kg
EBL 135	22.0cm	82cm	11.1kg
EBL 150	24.5cm	82cm	12.0kg
EBL 165	13.5cm	82cm	13.3kg
EBL 200	16.5cm	82cm	16.3kg
EBL 250	20.5cm	82cm	22.2kg
EBL 300	25.0cm	82cm	27.3kg
EBL 350	29.0cm	82cm	28.6kg

\* Recommended to check your local regulations and bylaws for any restrictions.

\*\* Recommended Max. Heights are based on assisted access with a gradient of; 1 in 6 for ramps up to 150cm. & 1 in 12 for ramps over 150cm.

\*\*\* Product weights are approximate.



Carry bag for Personal and Multipurpose models



Marked 'TOP' to show head of ramp, grip tape, pads at both ends



Multi-paneled, hinge system folds with a concertina action, allowing it to be folded to a single compact unit






Backrest width



Seat depth



Seat to Floor  
height

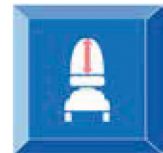

Backrest height



Driving unit width



Total length



Empty weight



Max user weight



Battery capacity



Max climbable  
obstacle height



Ground clearance



Max. safe slope



Motor capacity



Max. driving  
distance



Max. driving  
speed



The heaviest part

All information quoted is believed to be correct at time of print.  
Invacare® reserves the right to alter product specifications without prior consultation.



“ When you combine the power of **innovation** and **care**, the result is **Invacare.** ”



Unit 8, 110 - 120 Silverwater Road  
SILVERWATER NSW 2128  
Tel: 02 9649 2111  
E: [Sales@activemobility.com.au](mailto:Sales@activemobility.com.au)

[www.activemobility.com.au](http://www.activemobility.com.au)



**Active Mobility**