



Pressure Care



Active Mobility

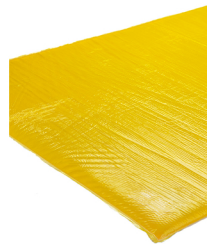

Pressure Bed Care

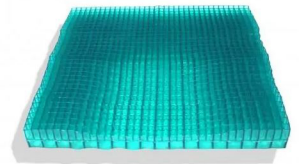

Pressure Seating Care




Personal Pressure Care

Action Mattress Overlay

Pressure Bed Overlay



- Made of 100% Akton polymer
- Used by hospitals and long-term aged care facilities to prevent and treat decubitus ulcers
- AKTON polymer will not migrate like gels, does not promote bacterial growth and is ideal for multiple client use
- Fluid proof
- Stain Resistant

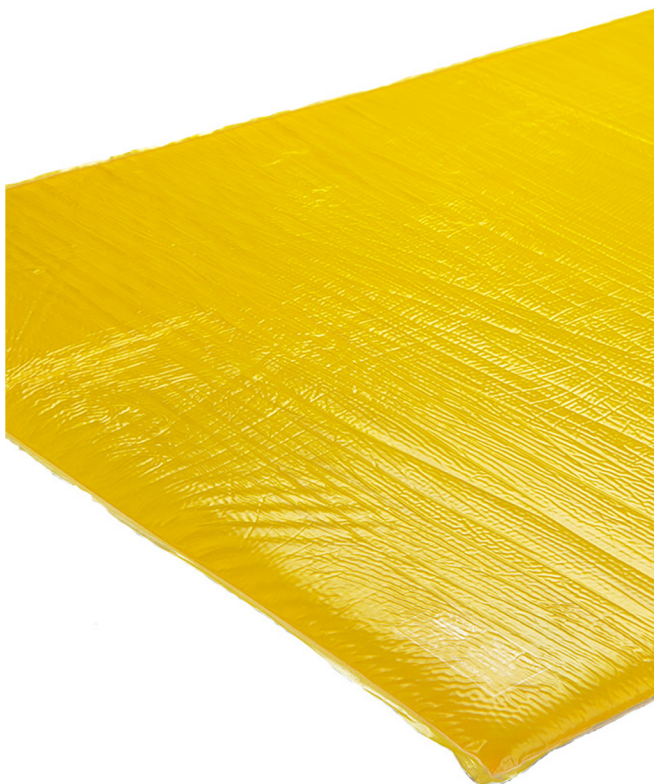


ACTION Overlay

Pressure Bed Overlay

- Quality AKTON® polymer mattress pad
- Promotes even pressure distribution
- AKTON® polymer gel, unlike fluid gel, does not displace from areas requiring pressure relief
- Isolates and protects areas where shear forces and pressure points are prevalent
- AKTON® polymer provides significant shock, vibration and impact absorption
- Pad is sealed by a waterproof film
- Pad does not promote bacterial growth, making it ideal for multi-client use

Specifications

Width	430 - 685mm
Length	658 - 2030mm
Height	22mm
Max User Weight	Unlimited

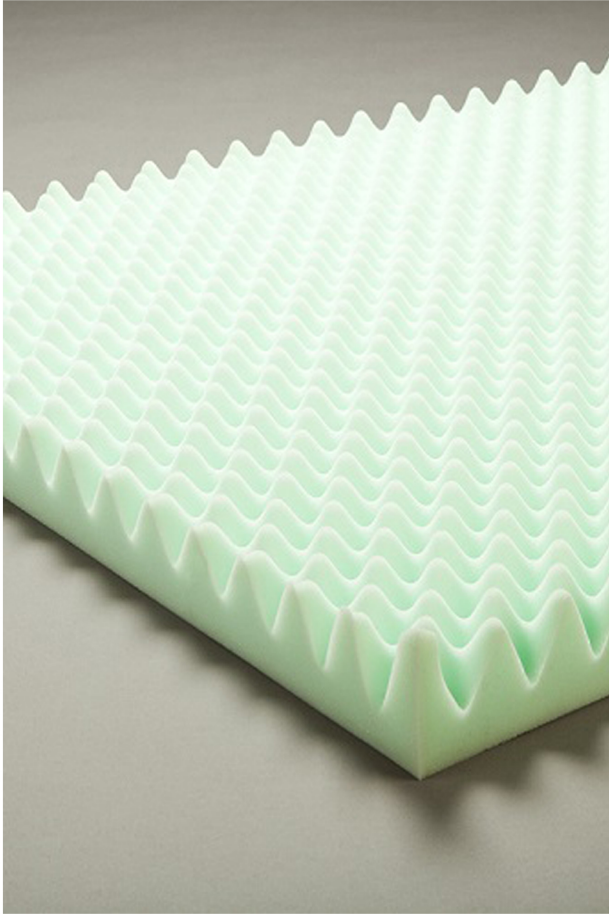

Convoluted Foam Overlay

Pressure Bed Overlay

- Soft open cell foam
- Provides additional comfort to existing mattress
- Convoluted profile provides increased air circulation

Specifications

Width	875 - 1485mm
Length	1850 - 1995mm
Height	75mm
Max User Weight	100kg



Opal Overlay

Pressure Bed Overlay

- Alternating pressure air overlay system
- Assists in the prevention of pressure problems for low risk patients
- Rows of air cell bubbles alternately inflate and deflate to assist in user comfort and to increase blood circulation
- Adjustable air pressure
- One piece PVC construction is hygienic and easy to clean
- Quiet operating pump hooks onto bed frame
- Single bed size

Specifications

Width	900mm
Length	2000mm
Height	65mm
Max User Weight	80kg



Prodigy Overlay

Pressure Bed Overlay

- Air mattress overlay system
- Ideal for clients who require low level pressure care
- Bed size is divided into three sections which are fastened together by press studs
- Polyvinyl cells are interconnected and are designed to help improve and maintain blood flow in order to prevent or treat skin issues
- Flexible cells adapt to patient contours and reduce the effects of shear forces
- Air pressure is adjustable and can be customised to accommodate client's individual needs
- Can be cleaned with soap and water
- Supplied with pump
- Single bed size

Specifications

Width (per section)	890mm
Length (per section)	680mm
Height	90mm
Max User Weight	135kg



ROHO Levelling Pad

Pressure Bed Overlay

- Quality foam and air mattress overlay system
- Open cell foam pad has an internal cut-out to accommodate a ROHO mattress section (supplied separately)
- Ideal where high level care needs are concentrated in a smaller area on the client
- Air cells are interconnected and their unique shape promote blood circulation to facilitate pressure prevention, healing and relief
- Air pressure is adjustable and can be customised to suit the user
- Durable neoprene construction is soft to suit skin sensitivity, and is easily cleaned
- Flexible cells adapt to patient contours and reduce shear forces while transferring or repositioning
- ROHO section is supplied with a pump

Specifications

Width	880mm
Length	1960mm
Height	80mm
ROHO Section Width	510mm
ROHO Section Length	860mm
Max User Weight	Unlimited

PRESSURE BED CARE




ROHO Mattress System *Pressure Bed Overlay*

- High quality air mattress overlay system
- Bed size is divided into four sections (ROHO or Foam) suit individual user requirements
- Mattress sections are easily fastened together by press studs
- Cells are interconnected and their unique shape promote blood circulation to facilitate pressure prevention, healing and relief
- Air pressure is adjustable and can be customised to suit the user
- Durable neoprene construction is soft to suit skin sensitivity, and is easily cleaned
- Flexible cells adapt to patient contours and reduce shear forces while transferring or repositioning
- Foam section, with it's waterproof cover, can be utilised where risk of pressure is not as critical
- ROHO sections are supplied with pump

Specifications

Width	860mm
Length (Per Section)	510mm
Height	90mm
Max User Weight	Unlimited

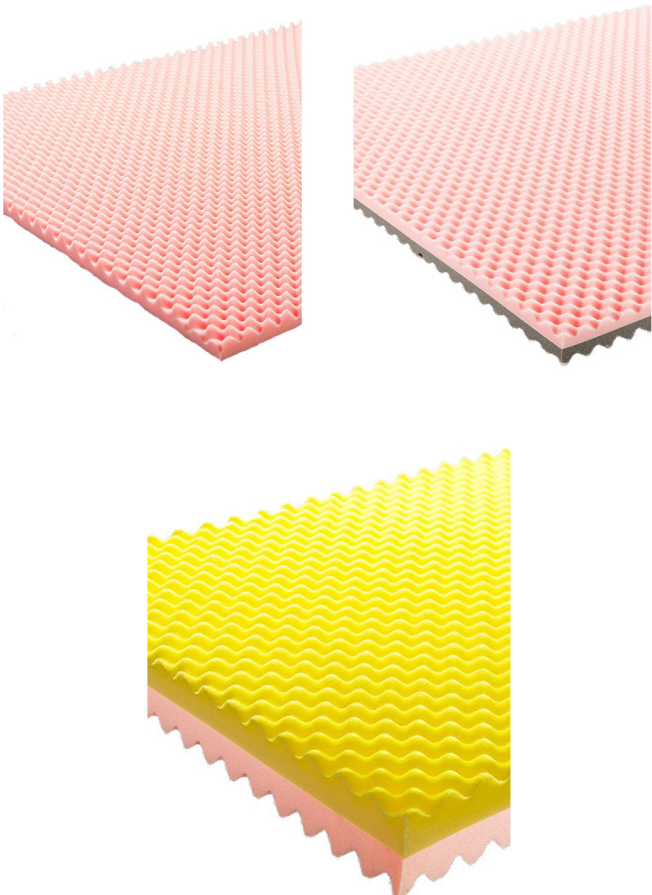

SAF Overlay *Pressure Bed Overlay*

- Gel permeated shock absorbent foam
- Open cell convoluted profile provides increased air circulation
- Memory foam contours to body shape to provide comfort and uniform pressure relief
- Anti-bacterial and machine washable up to 95°C
- Made in Denmark
- Single bed size

Specifications

Pink Model

Width	900mm
Length	2000mm
Height	40mm
Max User Weight	100kg

Pink/Grery Model

Width	900 - 1530mm
Length	2000mm
Height	40mm
Max User Weight	130kg

Pink/Yellow Model

Width	900mm
Length	2000mm
Height	60mm
Max User Weight	80kg



Silicone Fibre Overlay

Pressure Bed Overlay

- Silicone coated fibre filled overlay
- Cylindrical hollow core fibres for superior comfort and low level pressure relief
- Soft breathable cotton top side
- Water resistant polyurethane underside
- Elasticised corner straps hold overlay in position on mattress
- Machine washable

Specifications

Width	900mm
Length	1900mm
Height	100mm
Max User Weight	100kg



Sofflex Overlay

Pressure Bed Overlay

- Air mattress overlay system
- Ideal for clients who have moderate pressure care needs
- Bed size is divided into three sections which are fastened together by press studs
- Polyurethane cells are interconnected and are designed to help improve and maintain blood flow in order to prevent or treat skin issues
- Flexible cells adapt to patient contours and reduce the effects of shear forces
- Air pressure is adjustable and can be customised to accommodate client's individual needs
- Can be cleaned with soap and water
- Supplied with pump
- Single bed size

Specifications

Width (Per Section)	920mm
Length (Per Section)	700mm
Height	90mm
Max User Weight	135kg



Medical Sheepskin


Pressure Bed Overlay

- Genuine Australian lambswool overlay
- Provides additional comfort to existing mattress or chair
- Soft pile allows air circulation to keep skin dry
- Natural fibre reduces the risk of skin irritation
- Machine washable

Specifications

Width	570 - 640mm
Length	860 - 960mm

EquaGel General Cushion *Pressure Chair Cushion*

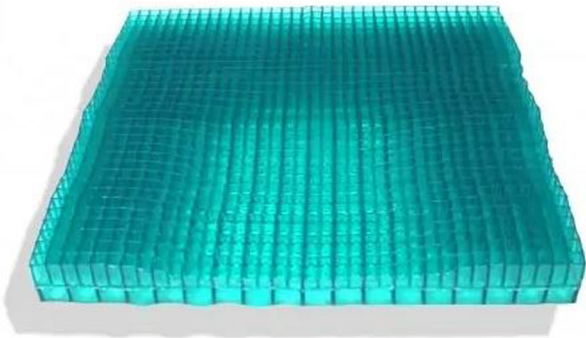

- 100% EquaGel
- 2" tall EquaGel in two stages
- The top layer is thin wall gel and a tight grid for a comfort layer
- The bottom layer is thicker wall gel and a wide grid for deep posture support
- The pelvic area of the cushion is contoured to more evenly distribute a user's weight than flat cushions, while retaining 1.25" of buckling gel under the pressure points
- Fluid resistant, vapor permeable and flame resistant cover

Specifications

SWL	220kg
Weight	3kg

EquaGel Protector Cushion *Pressure Chair Cushion*



- 100% EquaGel
- 2.5" tall EquaGel in two stages
- The top layer is a thin wall gel with a tight grid for a comfort layer
- The bottom layer has a thicker wall gell and a wide grid for deep posture support
- The pelvic area of the cushion is contoured to more evenly distribute a user's weight than flat cushions, while retaining 1.75" of buckling gel under the pressure points
- Fluid resistant, vapor permeable and flame resistant cover

Specifications

SWL	270kg
Weight	3kg

EquaGel Adjustable Protector Cushion

Pressure Chair Cushion




- 100% EquaGel
- 2.5" tall EquaGel in two stages
- The top layer is a thin wall gel with a tight grid for a comfort layer
- The bottom layer has thicker wall gel and a wider grid for deep posture support
- The pelvic area of the cushion is contoured to more evenly distribute a user's weight than flat cushions, while retaining 1.75" of buckling gel under the pressure points.
- Fluid resistant, vapor permeable and flame resistant cover
- Included are two different pelvic contour inserts that can increase the amount of therapeutic gel under the pelvis from 1.75" to 2", 2.25" or 2.5" deep.
- The gel contour inserts are slipped into the fabric pocket inside the cover, beneath the gel. You may customise your cushion by using no insert, only the thin, only the thick or both inserts,

Specifications

SWL	270kg
Weight	4kg

ACTION Adaptive Cube Pad

Pressure Chair Cushion

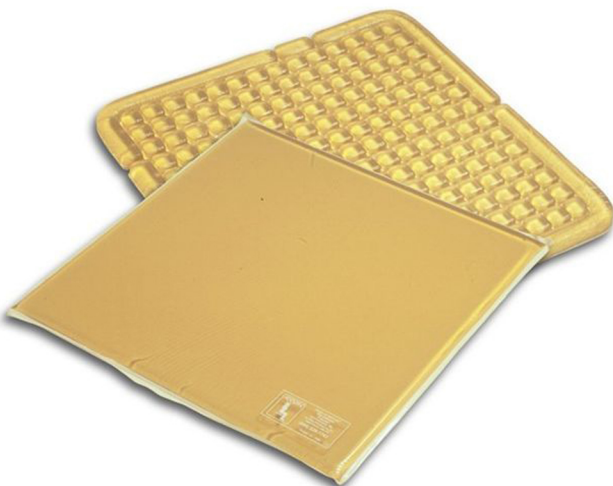

- Multifunctional applications for traditional and non-traditional seating challenges
- Designed to be cut and used for spot protection in product such as arm pads, foot plates, head/leg rests and splints/braces




Relax Easy Cushion

Pressure Chair Cushion

- Combination of layers of viscoelastic foam and FloGel
- FloGel promotes even pressure distribution while foam moulds to body shape providing additional comfort
- Cushion is fully rotational making it ideal for situations where low maintenance is required
- Cover is non-slip, fluid resistant and easy to clean

Specifications

Height	60mm
Max SWL	110kg

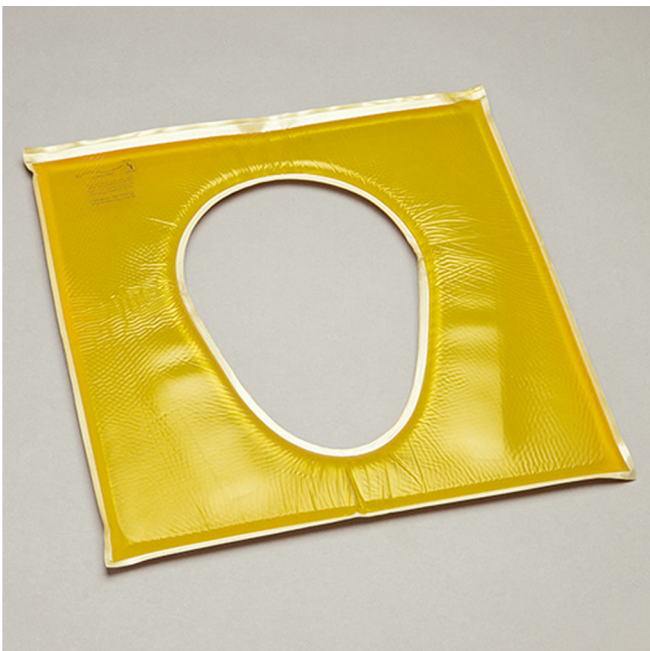

ACTION Polymer Commode Ring Cushion

Pressure Chair Cushion

- Quality AKTON® polymer
- Helps promote even pressure distribution
- AKTON® polymer, unlike fluid gel, does not displace from areas requiring pressure relief
- Isolates and protects areas where shear forces and pressure points are prevalent
- AKTON® polymer provides significant shock, vibration and impact absorption
- Pad does not promote bacterial growth, making it ideal for multi-client use
- Pad is sealed by a waterproof film
- Large model supplied as open front

Specifications

Length	41 - 46cm
Width	41 - 46cm
Height	16 - 25mm
Max User Weight	Unlimited



Silicone Fibre Wheelchair Pad

Pressure Chair Cushion

- Silicone coated fibre filled wheelchair cushion
- Manufactured from cylindrical hollow core fibres for greater comfort and shear reduction
- Back and arm pads provide greater user comfort
- Soft breathable cotton top side
- Fluid resistant PVC underside
- Machine washable

Specifications

Width	460mm
Depth	460mm
Height	75mm
Back Height	500mm
Arm Height	300mm
Max User Weight	100kg



Water Gel Foam Cushion

Pressure Chair Cushion


- Liquid gel filled cushion
- Water gel combination moulds around the contours of the user, assisting in even pressure distribution
- Gel is surrounded by foam which gives additional comfort and support
- Four way stretch vinyl on seating area reduces friction and shear
- Removable cover for cleaning
- Heavy duty vinyl sides reinforce the cushion
- Attached handle for easy transportation

Specifications

Width	440mm
Depth	440mm
Height	80mm
Max User Weight	110kg

Neptune Cushion

Pressure Chair Cushion

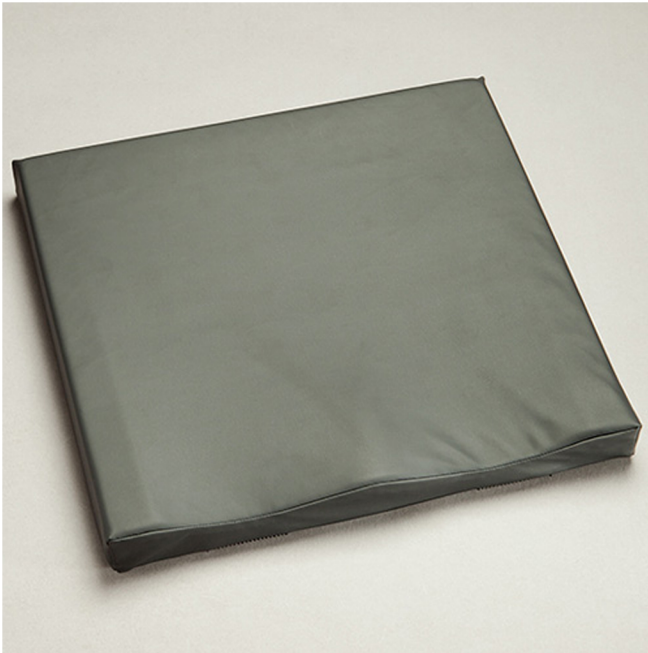

- Durable viscoelastic foam cushion
- Ideal for those with mild postural needs who are at low risk of skin breakdown
- Foam moulds to body shape, providing increased surface contact for improved pressure distribution and comfort
- Mildly contoured seating area provides increases postural positioning and stability
- Fluid resistant cover is zippered and features a non slip underside

Silicone Fibre Seat Pad

Pressure Chair Cushion



- Silicone coated fibre filled cushion
- Manufactured from cylindrical hollow core fibres for greater comfort and sear reduction
- Breathable cotton top side
- Fluid resistant PVC underside
- Machine washable

Specifications

Width	450mm
Depth	450mm
Height	75mm
Max User Weight	100kg

PRESSURE CHAIR CARE




Obus Forme Seat Support Pressure Chair Cushion

- Lightweight polyurethane foam cushion
- Ideal when sitting in a vehicle or when travelling
- Encourages proper alignment of the pelvis and thighs, evenly distributing body weight for extended sitting comfort
- Flexible sections to suit a wide range of seating surfaces
- Removable front section is easily unzipped to allow for use with narrow chairs and benches

Specifications

Width	450mm
Depth	420mm
Height	40mm
Max User Weight	100kg

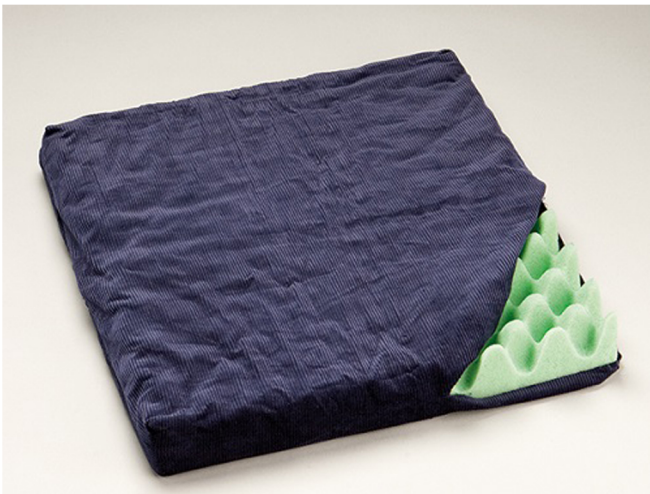

All Purpose Cushion Pressure Chair Cushion

- Soft open cell foam
- Convoluted surface for increased air circulation
- Corduroy fabric cover is zippered and removable for cleaning

Specifications

Width	460mm
Depth	460mm
Height	60mm
Max User Weight	100kg

PRESSURE CHAIR CARE



Coccyx Wedge

Pressure Care Wedge

- Soft open cell foam cushion with cutout area at rear that assists pressure relief on the tail bone (coccyx)
- Removable zippered cotton cover

Specifications

Width	390mm
Depth	400mm
Height	70mm
Max User Weight	100kg



Posture Wedge

Pressure Care Wedge

- Soft open cell foam cushion
- Designed to maintain a normal lumbar curve while sitting
- Promotes correct posture
- Anteriorly tilts pelvis to assist in straightening back
- Heavy duty zippered cotton cover

Specifications

Width	410mm
Depth	390mm
Height	30 - 75mm
Max User Weight	100kg



Sitt-on Wedge

Pressure Care Wedge

- Soft open cell foam cushion
- Designed to maintain a normal lumbar curve while sitting
- Promotes correct posture
- Anteriorly tilts pelvis to assist in straightening back
- Removable zippered cotton cover

Specifications

Width	400mm
Depth	330mm
Height	25 - 90mm
Max User Weight	100kg



Air Ring Cushion

Pressure Care Cushion

- Inflatable rubber cushion
- Ideal for short-term user requiring comfort and support
- Air pressure can be adjusted

Specifications

Outside Diameter	400mm
Inside Diameter	140mm
Height	80mm
Max User Weight	100kg



Convoluted Foam Ring Cushion

Pressure Care Cushion

- Soft open cell foam cushion
- Designed to maintain a normal lumbar curve while sitting
- Promotes correct posture
- Anteriorly tilts pelvis to assist in straightening back
- Removable zippered cotton cover

Specifications

Outside Diameter	430mm
Inside Diameter	150mm
Height	100mm
Max User Weight	100kg



Silicone Fibre Commode Ring Cushion

Pressure Care Cushion

- Silicone coated fibre filled cushion
- Manufactured from cylindrical hollow core fibres for greater comfort and shear reduction
- Made from fluid resistant PVC
- Adjustable velcro securing straps

Specifications

Width	410mm
Length	410mm
Height	75mm
Max User Weight	100kg



Therapod Back Support

Pressure Care Back Support

- Rigid framed orthopaedic back support
- Contoured support relieves back pain by realigning the spine and supporting correct posture
- Backrest contour can be altered via adjustable tension straps
- Covered in soft comfortable foam
- Adjustable nylon strap secures back support to chair or backrest
- High back model provides extended neck and head support for those who spend long hours driving or where vehicle back and headrest support is limited
- Blue brush-nylon removable cover

Specifications

Width	450mm
Height	600 - 770mm



Obus Forme Back Support

Pressure Care Back Support

- Rigid framed orthopaedic back support
- Contoured support relieves back pain by realigning the spine and supporting correct posture
- Covered in soft comfortable foam
- Supplied with removable adjustable lumbar pad
- Supplied with an elasticised strap to secure back support to chair or backrest
- High back model provides extended neck and head support for those who spend long hours driving or where vehicle back and headrest support is limited

Specifications

Width	450 - 500mm
Height	530 - 770mm



Apollo Back Support

Pressure Care Back Support

- Medium density foam back support cushion
- Provides comfort and support for the back and spine while sitting
- Maintains correct posture by supporting sideways movement
- Ideal for persons who sit for long periods
- Removable zippered cotton cover

Specifications

Width	360mm
Height	320 - 330mm



Swedish Back Support

Pressure Care Back Support

- Steel framed back support
- Provides comfort and support for the back and spine while sitting
- Maintains correct posture by supporting sideways movement
- Ideal for persons who sit for long periods
- Open mesh allows for increased air circulation
- Adjustable nylon buckled strap secures back support to chair or backrest

Specifications

Width	425mm
Height	390mm

ACTION Donut Head Pad *Pressure Care Headrest*



- Conform's to the patient's head for maximum support and pressure prevention
- Can be used in a variety of surgical procedures
- Dimensions: 20.5 x 7.5 x 4.5cm



ACTION Heel Protector Pad *Pressure Care Heel Pad*

- Should be used for any procedure that is lengthy and involves pressure on patient's heel area
- Easy to use
- Dimensions: 32 x 20 x 1.3cm

ACTION Heel Support

Pressure Care Heel Support

- Provides same protection as the heel protector pads, however, it is configured to secure and protect the heel area without the hook and loop attachment.
- Dimensions: 15 x 8.3 x 4.5cm



ACTION Horseshoe Head Pad

Pressure Care Heel Pad

- Low profile head pad used to protect the patient's ear while lying in a lateral position

Specifications

Outer Dimension	8.3 - 20.5cm
Inner Dimension	3.8 - 7.5cm
Height	2 - 7.5cm



PERSONAL PRESSURE CARE



ACTION Heel & Elbow Pads *Pressure Care Heel & Elbow Support*

- Recommended for patients 'at risk' or diagnosed with a compromised skin disorder
- Allows clients taking medications that affect skin integrity, or who are prone to skin tears or abrasions to continue daily activities or sleep comfortably



Sheepskin Slipper *Pressure Care Slippers*

- Provides comfort, warmth and assists in pressure care of the feet
- Genuine medical sheepskin
- Natural fibres reduces the risk of skin irritation
- Soft pile allows air circulation to keep dry
- Toe area is closed to increase comfort and support
- Adjustable velcro straps
- Soft sole allows slipper to be worn in bed
- Machine washable
- Supplied as pair

PERSONAL PRESSURE CARE



Heelift Suspension Boot

Pressure Care Boots

- Assists in the prevention and healing of pressure ulcers at the heel and ankle
- Delivers positioning and alignment of the lower limb
- Ideal for post-operative clients, for those who suffer from diabetes or persons who are at risk of developing foot drop
- An inbuilt bevelled foam wedge provides support for the calf, and ensures that the heel is suspended to reduce pressure
- Secured in place with soft adjustable Velcro straps
- Ventilation holes allow added air circulation
- Stiffened base maintains the boot's structure and supports the forefoot
- Low friction surface allows the boot to slide easily over the mattress surface, improving client mobility
- Supplied with an extra foam pad to provide added elevation
- Optional convoluted foam surface for sensitive skin
- Foam can be cut to suit individual requirements
- Latex-free medical grade foam is autoclavable machine, washable and dryer safe
- One size fits all
- Supplied as each



Silicone Fibre Heel Protector

Pressure Care Heel Protector

- Padded heel cushion
- Provides comfort and assists in the prevention of pressure build-up in the heel area
- Manufactured from millions of silicone coated cylindrical hollow cored fibres
- Silicone coating acts as a lubricant to eliminate friction and reduce shear forces
- Adjustable velcro securing straps
- Internal seams are carefully positioned to avoid creating pressure points
- Breathable cotton cover
- Machine washable
- One size fits all
- Supplied as each



Silicone Fibre Leg Protector

Pressure Care Leg Protector

- Padded leg cushion
- Provides comfort and assists in the prevention of pressure build-up in the leg area
- Manufactured from millions of silicone coated cylindrical hollow cored fibres
- Silicone coating acts as a lubricant to eliminate friction and reduce shear forces
- Adjustable velcro securing straps
- Internal seams are carefully positioned to avoid creating pressure points
- Breathable cotton cover
- Machine washable
- One size fits all
- Supplied as each



Peanut Pillow

Pressure Care Neck Pillow

- Soft contoured neck pillow that supports the head and neck
- Can be used either in bed or chair

Specifications

Width	300mm
Height	100mm

PERSONAL PRESSURE CARE



Push Cush

Award-winning inflatable foot cushion

- Award winning inflatable foot cushion and exercise aid
- Relieves swollen legs and feet when travelling or sitting for long periods
- Ideal for those suffering from fluid and circulatory problems
- Pedal/treadle action exercises promote circulation and oxygen to the legs and feet
- Works important muscles and keeps joints mobile
- Widely used to prevent DVT (Deep Vein Thrombosis) for those travelling in confined positions

Specifications

Width	400mm
Depth	330mm



SAF Heel Protector

Pressure Heel Protector

- Provides comfort and assists in the prevention of pressure build-up in the heel area
- Gel permeated shock absorbent foam
- Open cell convoluted profile provides increased air circulation
- Top soft layer of yellow foam is laminated with medium density grey foam
- Yellow foam moulds to body shape to provide optimum comfort while grey foam provides uniform pressure distribution and support
- Secured in place with soft adjustable Velcro straps
- Low friction surface allows the heel protector to slide easily over the mattress surface
- Anti-bacterial and machine washable
- Made in Denmark

PERSONAL PRESSURE CARE



Sheepskin Elbow Protector

Pressure Elbow Protector

- Provides comfort and assists in the prevention of pressure build-up in the elbow area
- Genuine medical sheepskin
- Soft pile allows air circulation to keep skin dry
- Natural fibres reduces the risk of skin irritation
- Adjustable velcro securing straps
- Machine washable
- One size fits all
- Supplied as a pair



Sheepskin Heel Protector

Pressure Heel Protector

- Provides comfort and assists in the prevention of pressure build-up in the heel area
- Genuine medical sheepskin
- Soft pile allows air circulation to keep skin dry
- Natural fibres reduces the risk of skin irritation
- Adjustable velcro securing straps
- Machine washable
- One size fits all
- Supplied as a pair



Silicone Fibre Elbow Protector

Pressure Heel Protector

- Padded elbow cushion
- Provides comfort and assists in the prevention of pressure build-up in the elbow area
- Manufactured from millions of silicone coated cylindrical hollow cored fibres
- Silicone coating acts as a lubricant to eliminate friction and reduce shear forces
- Adjustable velcro securing straps
- Internal seams are carefully positioned to avoid creating pressure points
- Breathable cotton cover
- Machine washable
- One size fits all
- Supplied as each