

Shower & Bath

Baria



Height adjustable bariatric electric shower trolley

SKU 4435028

The Baria shower trolley was designed to provide personal hygiene care to people with morbid obesity. The Baria shower trolley is the only one in the world to have a fifth motorized wheel to facilitate smooth movement in tight spaces, allowing it to be moved by a single caregiver. Comfort, stability and security are guaranteed to patients thanks to its studied center of gravity and its double high barriers in stainless steel. Its high chassis allows the Baria to be used with all transfer systems on the market.



Baria 3500.00 model fitted with wheel 360 ° electric «effortless».



Active Mobility

Shower & Bath

Baria

UNIQUE FEATURES



Ø150mm wheels



Emergency stop device



10 functions hand control



Sensitive push handle




Automatic locking



Removable battery

- Multiple positions to facilitate care.
- Bath position.
- 10cm stretching position.
- Fifth motorized wheel, 360 ° swivel to facilitate movement.
- Double high safety barriers in stainless steel.
- Wireless remote control.
- Lifting amplitude.
- Unbreakable tray.
- Capacity of 350kg.
- Ultra-light aluminum structure.
- Bi-centralized braking.

Height adjustable bariatric electric shower trolley

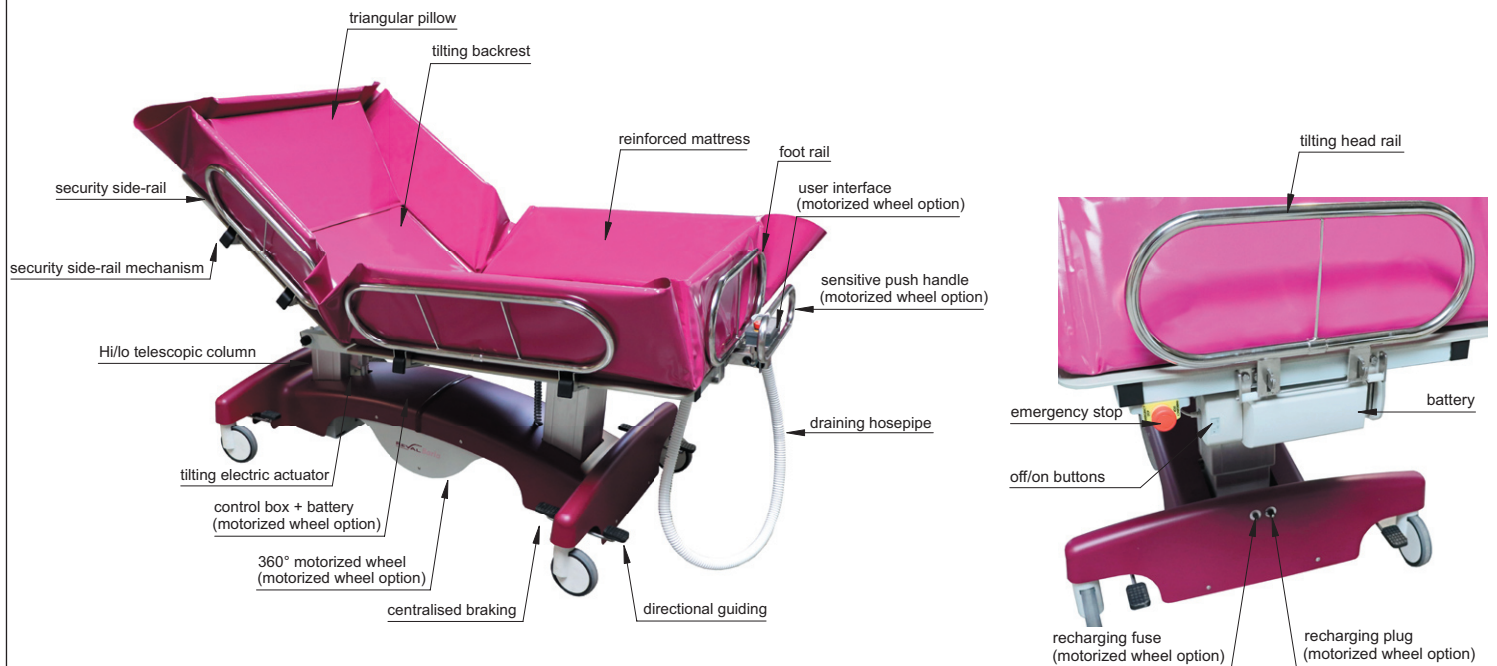

Shower & Bath

Baria



Wireless remote control option


360° electric wheel option

Panel 3011.01

Panel 3012.01

Height adjustable bariatric electric shower trolley

REF. 3500.00

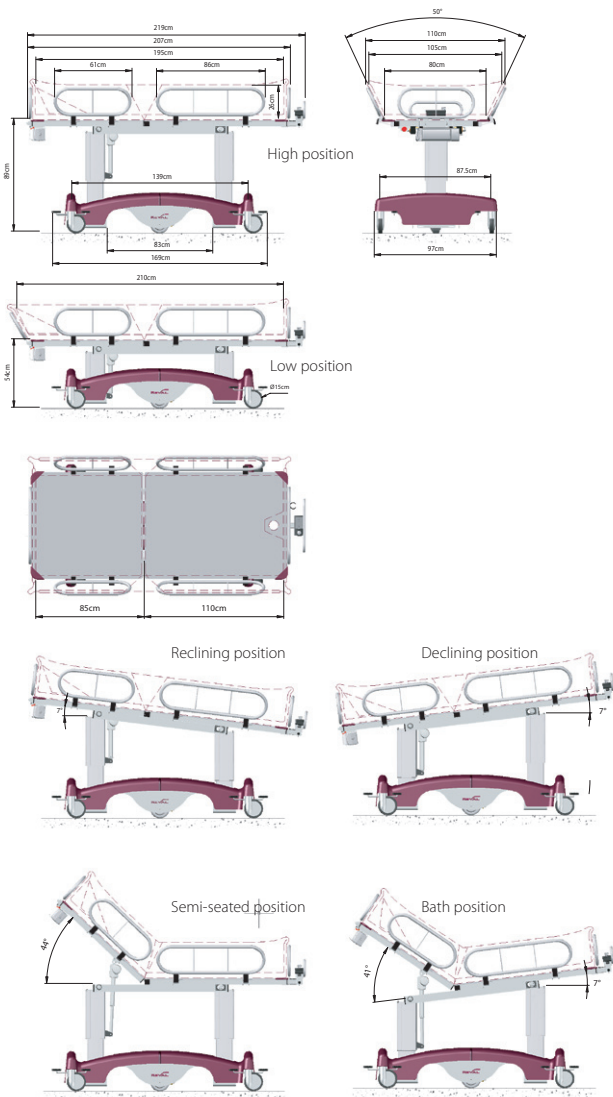

SPECIFICATIONS:

Chassis

- Aluminium & paint finished rolling chassis.
- ABS plastic chassis covers.
- x 4 Guiding and braking castors.
- Steering handle.

Platform

- x2 Articulating patient support sections in polypropylene.
- Stainless steel articulating hinges and fixations.
- x6 Stainless steel security side rails with locking functions.
- Tilting stainless steel head rail.

Mattress

- Soft PVC waterproof mattress designed for heavy usage, double thickness.
- Drainage outlet and waste water hose at foot end Ø40mm.
- Triangular headrest, in soft durable waterproof PVC material.

Electric components

- 24Vdc IPX6 rated variable height platform columns.
- 24Vdc IPX6 rated backrest actuator
- IPX6 Control module
- IPX6 10 function hand set with position memorization and charge indicators
- 5Ah / 24V Battery
- 24V battery and charger module.
- Emergency stop button.

Performances:

- SWL – Max load 350Kg
- Weight of appliance – 119Kg
- Weight of appliance including the electric wheel option – 141Kg
- Usage capacity - max 18 cycles @ 350Kg
- Usage capacity - max 26 cycles @ 260Kg
- Up/down speed 9mm/s @ 350kg
- Action radius of the electric wheel on full battery : 10 to 15km
- Max assisted speed with the electric wheel : 4.5km/h

Options & Accessories:

- 4435028-10 10-button wireless remote + receiver + support
- 4435028-20 «Effortless» 360° electric wheel assistance
- 4435028-30 Additional battery
- 4435028-40 Extra charger and base for batteries

Conformity:

- 93/42/CEE (Medical Equipment Directive)
- NF EN 12182 (Tech. support systems for less able persons)
- ISO 14971 (Medical Devices - Application of Risk)

