

IC777 PRODUCT PORTFOLIO







ICARE IC777

CARING
FOR YOUR
COMFORT



Active Mobility

IC777 PREMIUM HOMECARE BED



FABRIC


Heavy commercial grade velvet feel fabric.



Mist




Royal

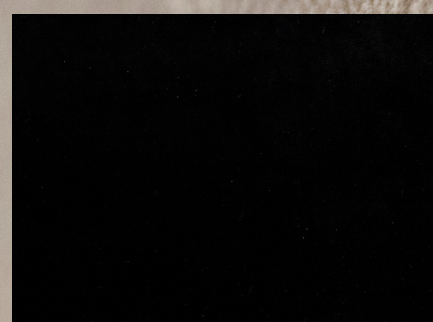


Noir

VINYL

Heavy commercial grade jacquard, vinyl fabric.

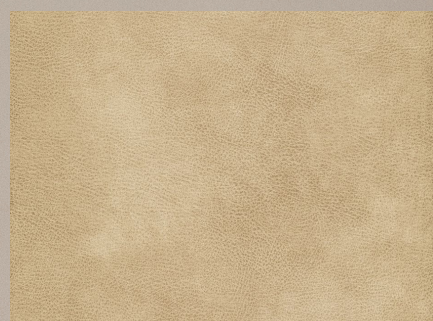

Sand

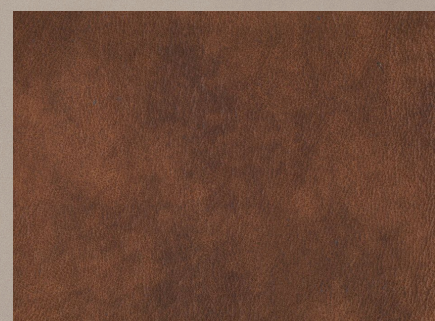

Bison



Slate

This four function bed, features dual lift motors totaling 12,000 newtons lifting power, unlimited under-bed clearance and eight heavy duty castors.



62° Head Lift



42° Knee Break



Height Adjustable
Min 22 - Max 66cm



Trendelenburg and
Reverse Trendelenburg

Bring beauty and character into your room with exceptional bed functions, made in Australia.

AMSTERDAM

Headboard and footboard



Pictured in Mist Fabric



COLOGNE

Headboard and footboard



Pictured in Mist Fabric



KARLSTEAD

Headboard and footboard



Pictured in Royal Fabric



IC700 PARTNER BED

By having a co-sleeping partner bed configuration this can significantly help with recovery and providing positive care. IC700 is a static bed, designed to sit right alongside the IC777 and looks just as elegant.

Available in Long Single or King Single sizes.



ACCESSORIES

Complete care at home! The IC777 is compatible with full range of accessories to make independence easier every day.



AFTER DARK FEATURES



Flash light on the back of the remote



Under bed LED lighting



Two personal programmable memory position settings



Zero gravity pre-set functions



Patient transfer pre-set function



Backlit remote

REMOTE

Backlit cordless remote.



- Backlit remote activated by a touch of a button.
- Remote preprogrammed functions:
 - M1 - user personal settings.
 - M2 - user personal settings.
 - ZG Zero Gravity - preset function-override for personal settings.
 - PT Patient Transfer - preset function-override for personal settings.
- LED under bed lighting.
- Flashlight on/off button.
- Indicator light: red- battery low, green-charged

BED CODES & SIZES

Sizes for every room size and lifestyle.

BED SIZE	CODE	MATTRESS SURFACE SIZE	BASE SIZE
Long Single	1481718	2030 x 900mm	2150 x 960mm
King Single	1481718	2030 x 1070mm	2150 x 1120mm
Long Double • One hand control	1481718	2030 x 1350mm	2150 x 1400mm
One Piece Queen • One hand control	1481718	2030 x 1500mm	2150 x 1580mm
Split Queen		2030 x 750mm (per side)	2150 x 1580mm
Split King		2030 x 900mm (per side)	2150 x 940mm • as partner

NOTES:

SPLIT QUEEN

- Separate head and knee functions
- One piece hi-lo and trendelenberg
- Individual hand control per side
- Fits a one piece queen headboard

SPLIT KING

Always remember that TWO Long Single beds together make a King size bed.

TECHNICAL SPECS

Key details and information.

- 200kg safe working load on long single and king single sizes and 300kg for double and queen sizes. The safe working load includes any model of Icare mattresses.
- Long-life rechargeable battery back up, lasting 2-5 days of usual use.
- Easy press button remote, backlit when activated for use. Wipeable.
- Hard wearing commercial grade fabric. Vinyl is treated with antibacterial properties.
- Strong metal mesh mattress platform allowing increased airflow
- Actuators - 1x 6000 newton - 24 watt head adjust actuator. 1x 6000 newton - 24 watt knee adjust actuator. 2x 6000 newton - 24 watt height adjust actuators (2x 8000 newton on double and queen sizes)
- No pull motors for added safety.
- Wheels - 8x heavy duty swivel castors with one lockable castor on each corner.



HEADBOARDS & FOOTBOARDS

Codes and measurements.



AMSTERDAM HEADBOARD

LONG SINGLE

Code: HB7AS

Width: 960mm

KING SINGLE

Code: HB7AKS

Width: 1120mm

DOUBLE

Code: HB7AD

Width: 1400mm

QUEEN

Code: HB7AQ

Width: 1580mm

AMSTERDAM FOOTBOARD

LONG SINGLE

Code: FB7AS

Width: 960mm

KING SINGLE

Code: FB7AKS

Width: 1120mm

DOUBLE

Code: FB7AD

Width: 1400mm

QUEEN

Code: FB7AQ

Width: 1580mm



COLOGNE HEADBOARD

LONG SINGLE

Code: HB7CS

Width: 960mm

KING SINGLE

Code: HB7CKS

Width: 1120mm

DOUBLE

Code: HB7CD

Width: 1400mm

QUEEN

Code: HB7CQ

Width: 1580mm

COLOGNE FOOTBOARD

LONG SINGLE

Code: FB7CS

Width: 960mm

KING SINGLE

Code: FB7CKS

Width: 1120mm

DOUBLE

Code: FB7CD

Width: 1400mm

QUEEN

Code: FB7CQ

Width: 1580mm



KARLSTEAD HEADBOARD

LONG SINGLE

Code: HB7KS

Width: 960mm

KING SINGLE

Code: HB7KKS

Width: 1120mm

DOUBLE

Code: HB7KD

Width: 1400mm

QUEEN

Code: HB7KQ

Width: 1580mm

KARLSTEAD FOOTBOARD

LONG SINGLE

Code: FB7KS

Width: 960mm

KING SINGLE

Code: FB7KKS

Width: 1120mm

DOUBLE

Code: FB7KD

Width: 1400mm

QUEEN

Code: FB7KQ

Width: 1580mm

ICARE IC777

A NEW
WAY TO
CARE



Active Mobility

IC777 PREMIUM HOMECARE BED







Active Mobility

ACTIVE MOBILITY SYSTEMS PTY LTD

Unit 8/110-120 Silverwater Road

Silverwater NSW 2128

sales@activemobility.com.au

www.activemobility.com.au