



Luna Ceiling Hoist

Fixed / Portable

The lightweight and versatile Luna ceiling hoist installs in seconds and when paired together with the compatible tracking systems, offers optimal room coverage.



Fixed

SKU:

2333997 - 200kg
2333995 - 275kg

Portable

SKU:

2333998 - 200kg
2333996 - 275kg

Features



Emergency Lower
Emergency lowering functions



Weight Capacity
200kg
275kg



Wet Environments
IPx4 tested for wet environments



Easy Charging
Charges through the hand control via the wall mounted docking station

Full Features

- + Has the ability to be used as either a fixed or portable ceiling hoist in one and the same unit.
- + Changing Luna tracks with a Luna portable can be done within seconds by simply detaching and reattaching the unit to the track.
- + Cost effective for institutions who need only invest in only one modular model to cover all their lifting requirements.
- + Emergency lowering functions on hand control & Luna motor.
- + Charging through the ergonomic hand control on wall mounted docking station.

Specifications

- + Max Capacity : 200/275 kg
- + Length of Lifting Strap : 2200mm
- + Weight : 7 kg (w/o spreader bar)
- + Battery Operation : 36 V lead acid batteries
- + Charging : 100-240 V, 5A
- + Chassis: Aluminium
- + Housing : High impact moulded ABS
- + Charging Time : Aprpr 4.5-5 hours to charge an empty battery
- + Lifting Speed Down : 44mm/sec (w 89kg load)
- + Lifting Speed Up : 30mm/sec (w 89kg load)
- + No. of Lifts per Charge : 50 lifts of 0.5m (w 100kg)
- + Protection Class : IPx4

Accessories


Hand Control IPx4
1320mm:
233265-1320
1620mm:
233265-1620
1900mm:
233265-1900

Hand Control IPx7:
2333280

3,4 Point spreader bars:
3 Point: 2333030
4 Point: 2333028

Weigh Scales:
233188
2333260

Extension Strap in various lengths also available.

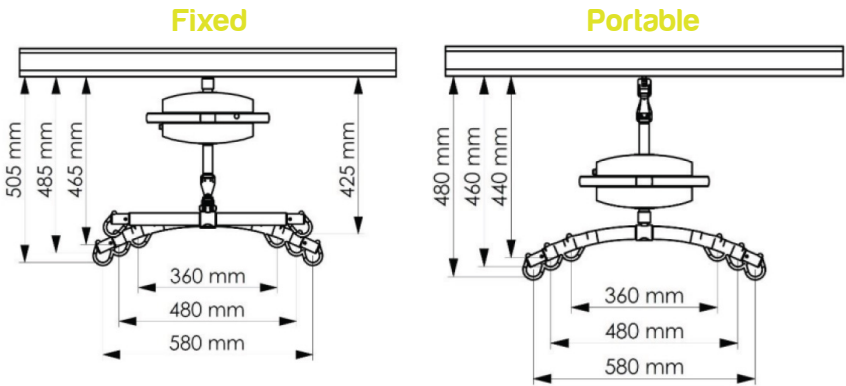

Unit 8, 110-120 Silverwater Rd
Silverwater NSW 2128
Ph: 02 9649 2111

activemobility.com.au

IMPORTANT : Always be sure to read your owner's manual and installation guide thoroughly before operation.

Warranty Information:
Hoist: 2 years
Batteries: 6 months

