

# PRIDE® Pathrider 140 XL

An exciting blend of power, precision and exceptional style, the Pathrider 140 XL is the perfect outdoor scooter.




## Colours Available



Yellow



Candy Apple Red

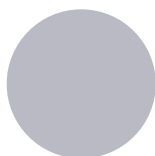

Silver

## Features



**Seat**  
High back captain seat with pillow top



**Wheels**  
4 Wheel



**Weight Capacity**  
181kg  
Bariatric suitable



**Lighting**  
LED headlights, front shield, indicators, hazard lights

# Full Features

- Powerful drivetrain with 24-Volt DC motor
- Full suspension for optimum handling
- Sealed hydraulic brakes
- Large 33cm pneumatic tyres for excellent outdoor performance
- Easy drive tiller with wraparound handles
- Easily accessible tie down points (for transport of unoccupied scooter)
- Rearview mirror
- Front shield
- High back pillow top seat with headrest
- Front & rear suspension
- Front basket
- LED lower headlight
- Hazard lights / indicators
- Sleek, sporty style

# Specifications

SPECIFICATIONS	PRIDE PATHRIDER 140XL
Maximum Weight Capacity	181kg
Maximum Speed <sup>1</sup>	15km/h
Maximum Range <sup>1,4</sup>	Up to 35km with 75AH batteries
Length <sup>2</sup>	142.8cm
Width <sup>2</sup>	69.2cm
Scooter Weight (without battery) <sup>2,5</sup>	117.70kg
Ground Clearance <sup>2</sup>	13.33cm
Seat to Ground	72.39cm - 68.58cm
Turning Radius <sup>2</sup>	Right: 210.8cm Left: 208.3cm
Battery Requirements <sup>3</sup>	2 x 12-Volt, Deep Cycle 75AH Batteries
Charger	8Amp, Off Board
Seating	Pillow Top, 21.77kg, Black Vinyl, Width 49.5cm, Depth 44.5cm
Motor	Rear Wheel Drive, Sealed Transaxle, 24-volt DC Motor

<sup>1</sup>Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tyre condition. These specifications can be subject to a variance of (+/- 10%).

<sup>2</sup>Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or - 3%).

<sup>3</sup>AGM or gel cell type required.

<sup>4</sup>Tested in accordance with ANSI/RESNA, WC Vol2, section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance.

Test conducted at maximum weight capacity.

<sup>5</sup>Battery weight may vary based on manufacturer.

## Warranty Information:

Scooters: 12 Month Warranty for all electronic parts incl controllers & battery charges

Scooter Frames: 24 Month warranty on the frame incl platform, fork, seat post, drivetrain, differential, motor, brakes and welds.

Batteries: 6 Month warranty

**IMPORTANT :** Always be sure to read your owner's manual thoroughly before operating your machine.

**WARNING :** Do not exceed 10km/ph on public roads or paths.

### Silverwater

Unit 8, 110-120 Silverwater Rd  
Silverwater NSW 2128  
Ph: 9649 2111

### Newcastle

7 Poynton Pl  
Thornton, NSW 2322  
Ph: 4047 3608